

General

Name of Run	230719_CombinedStd_LessTimeR3.run	Name of Analysis	NOx_P04.ANL
Date of Report	7/19/2023	System	AA500
Run Start	7/19/2023 4:39:07 PM	Run Stop	6:38:55 PM
Operator Run/Rec.	CLMR / CLMR	Software Version	AACE 8.02 SP2
Comment	Recalculate from Run230719_CombinedStd_LessTime.run		

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Slope only	Slope only	
Corr.Coeff.(r)	0.9995	0.9997	
AUFS/Gain	0.45/22	0.10/99	Data
Offset	-15845	3635	Field
Sensitivity	0.3718	0.0852	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	0.000	0.000	0
1	921	P	Primer	0.191	0.203	1
2	921	D	Drift	0.192	0.203	1
3	921	C	0.200	0.199	0.204	1
4	901	C	0.100	0.103	0.106	1
5	902	C	0.050	0.052	0.059	1
6	903	C	0.025	0.027	0.030	1
7	904	C	0.010	0.011	0.014	1
8	905	C	0.005	0.006	0.010	1
9	906	C	Blank (Ac. Seawater)	-0.003	0.004	1
10	0	B	Baseline	0.000	0.000	1
11	921	H1	High	0.271	0.204	1
12	905	L1	Low	0.005	0.010	1
13	905	L1	Low	0.006	0.010	1
14	906	S	Blank	-0.005	0.004	1
15	907	K	RecoveryBlank spiked	0.002	0.007	1
16	908	QC1	CRM N03	0.018 F	0.008	1
17	909	QC2	CRM P04	-0.006	0.010	1
18	1	S	37	0.001	0.003	1
19	2	S	N03_5	0.001	0.001	1
20	3	S	N03_6	0.006	0.001	1
21	4	S	P04_3	-0.004	0.001	1
22	5	S	P04_4	-0.004	0.006	1
23	6	S	P04_5	0.000	0.001	1
24	7	S	43	0.003	0.003	1
25	8	S	23	0.000	0.006	1
26	9	S	32	-0.011	0.002	1
27	10	S	29	0.000	0.007	1
28	11	K	29 spiked	0.000	0.005	1
29	908	QC1	CRM N03	0.031 F	0.008	1
30	909	QC2	CRM P04	0.000	0.011	1
31	0	B	Baseline	0.000	0.000	1
32	921	D	Drift	0.192	0.203	1
33	12	S	21	0.000	0.002	1
34	13	S	31	0.000	0.002	1
35	14	S	27	0.000	0.001	1
36	15	S	35	0.000	0.002	1
37	16	S	36	0.000	0.004	1
38	17	S	33	0.000	0.001	1
39	18	S	34	0.006	0.002	1
40	19	S	40	0.010	0.003	1
41	20	S	52	0.000	0.003	1
42	21	S	N03_3	-0.010	0.002	1

43	22	K	N03_3 spiked	0.050	0.004	1
44	908	QC1	CRM N03	0.000 F	0.008	1
45	909	QC2	CRM P04	-0.005	0.010	1
46	0	B	Baseline	0.000	0.000	1
47	921	D	Drift	0.192	0.203	1
48	23	S	38	0.009	0.005	1
49	24	S	39	0.016	0.001	1
50	25	S	44	0.000	0.006	1
51	26	S	45	-0.007	0.001	1
52	27	S	41	-0.009	0.001	1
53	28	S	47	-0.002	0.002	1
54	29	S	25	-0.009	0.001	1
55	30	S	49	0.038	0.001	1
56	31	S	46	0.010	0.002	1
57	32	S	42	0.005	0.004	1
58	33	K	42 spiked	0.056	0.007	1
59	908	QC1	CRM N03	-0.007 F	0.008	1
60	909	QC2	CRM P04	-0.004	0.010	1
61	0	B	Baseline	0.000	0.000	1
62	921	D	Drift	0.192	0.203	1
63	34	S	50	0.002	0.003	1
64	35	S	51	0.000	0.003	1
65	36	S	55	0.000	0.004	1
66	37	S	68	0.000	0.001	1
67	38	S	59	-0.008	0.006	1
68	39	S	65	0.000	0.004	1
69	40	S	67	0.000	0.004	1
70	41	S	58	0.000	0.005	1
71	42	S	54	0.003	0.002	1
72	43	S	64	0.051	0.004	1
73	44	K	64 spiked	0.042	0.026	1
74	908	QC1	CRM N03	0.000 F	0.007	1
75	909	QC2	CRM 904	0.000	0.009	1
76	0	B	Baseline	0.000	0.000	1
77	921	D	Drift	0.192	0.203	1
78	45	S	61	-0.007	0.005	1
79	46	S	62	0.000	0.006	1
80	47	S	56	-0.003	0.003	1
81	48	S	66	0.000	0.005	1
82	49	S	48	-0.005	0.002	1
83	50	S	57	0.001	0.002	1
84	51	S	53	0.021	0.004	1
85	52	K	53 spiked	0.049	0.013	1
86	908	QC1	CRM N03	0.000 F	0.007	1
87	909	QC1	CRM P04	0.000 F	0.009	1
88	0	B	Baseline	0.000	0.000	1
89	921	D	Drift	0.192	0.203	1
90	0	B	Final Base	0.000	0.000	1

Corrections

Channel	1	3
Baseline	done	done
Drift	done	done
Carryover	done	done
%	0.62	0.22

Spike Recovery

RecoveryBlank spiked Added	0.002	0.002
RecoveryBlank spiked Found	0.007	0.003
Recovery %	271.248	180.433
29 spiked Added	0.002	0.002
29 spiked Found	0.000	-0.001
Recovery %	0.000	-95.345
N03_3 spiked Added	0.003	0.002
N03_3 spiked Found	0.060	0.002
Recovery %	2058.321	117.579
42 spiked Added	0.002	0.002
42 spiked Found	0.050	0.003
Recovery %	3069.446	163.961
64 spiked Added	-0.002	0.002
64 spiked Found	-0.009	0.022
Recovery %	405.357	1270.019


53 spiked Added	0.000	0.002
53 spiked Found	0.028	0.009
Recovery %	8329.823	518.520

Calibration

Channel 1 NOx

Type	Observed	Calculated	Target	% Error
1C	59041	0.199	0.200	-0.72
2C	32327	0.103	0.100	3.29
3C	18023	0.052	0.050	4.55
4C	10829	0.027	0.025	6.47
5C	6358	0.011	0.010	6.71
6C	4980	0.006	0.005	15.13
7C	2397	-0.003	0.000	-


Calibration Fit: Slope only
 Correlation Coefficient: 0.9995
 Carryover (%): 0.6217
 Calibration Equation: $y = bx$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 0.000000E+000
 b = 2.803822E+005

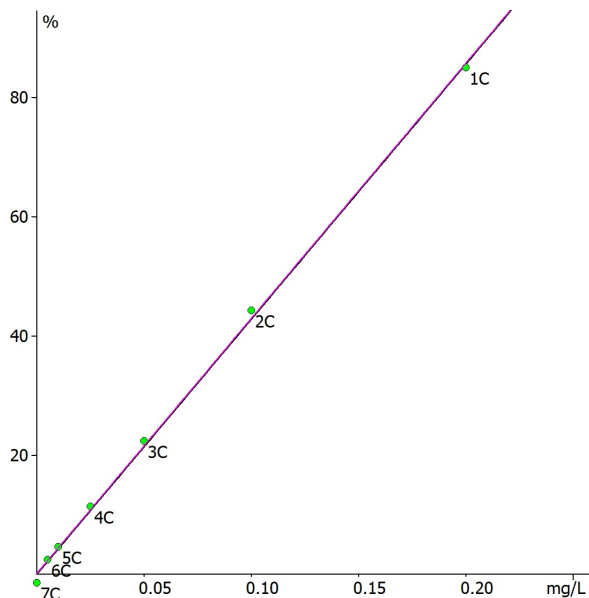

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58744	0.204	0.200	1.99
2C	32002	0.106	0.100	5.73
3C	19217	0.059	0.050	17.50
4C	11339	0.030	0.025	19.22
5C	7061	0.014	0.010	40.87
6C	5902	0.010	0.005	96.57
7C	4200	0.004	0.000	-

Calibration Fit: Slope only
 Correlation Coefficient: 0.9997
 Carryover (%): 0.2194
 Calibration Equation: $y = bx$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 0.000000E+000
 b = 2.721630E+005

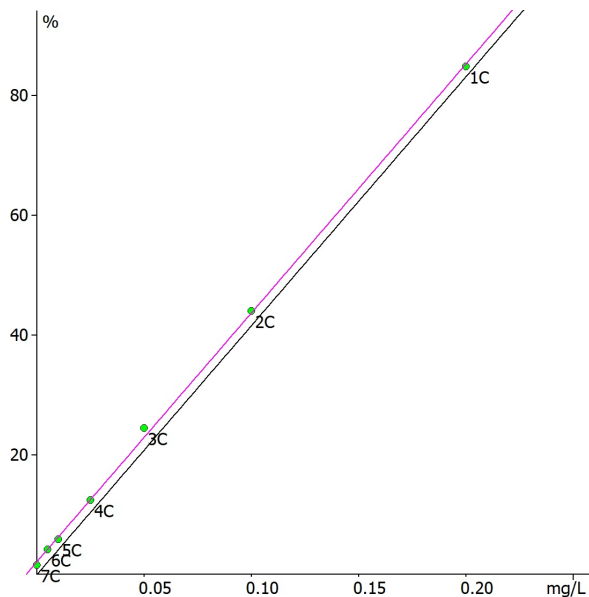

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero