

General

Name of Run	231109_Q2_NP2.run	Name of Analysis	NOx_P04.ANL
Date of Report	11/9/2023	System	AA500
Run Start	11/9/2023 1:06:45 PM	Run Stop	2:49:28 PM
Operator Run/Rec. CLMR		Software Version	AACE 8.02 SP2
Comment			

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	1.0000	1.0000	
AUFS/Gain	0.42/24	0.10/99	Data
Offset	55	5055	Field
Sensitivity	0.3499	0.0841	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.002	-0.002	0
1	921	P	Primer	0.201	0.199	1
2	921	D	Drift	0.200	0.200	1
3	921	C	0.200	0.200	0.200	1
4	901	C	0.100	0.101	0.101	1
5	902	C	0.050	0.050	0.049	1
6	903	C	0.025	0.025	0.024	1
7	904	C	0.010	0.010	0.011	1
8	905	C	0.005	0.004	0.005	1
9	906	C	Blank	0.000	0.000	1
10	0	B	Baseline	-0.002	-0.002	1
11	921	H1	High	0.199	0.199	1
12	905	L1	Low	0.004	0.005	1
13	905	L1	Low	0.004	0.005	1
14	906	S	Blank	0.000	0.000	1
15	907	K	Blank spiked	0.138	0.137	1
16	908	QC1	CRM N03 [0.030] Old	0.010 F	-0.001	1
17	909	QC2	CRM N03 [0.030] New	0.029	0.001	1
18	1	S	243	0.004	0.002	1
19	2	S	254	0.006	0.002	1
20	3	S	265	0.007	0.003	1
21	4	S	276	0.007	0.002	1
22	5	S	287	0.012	0.002	1
23	6	S	298	0.014	0.003	1
24	7	S	309	0.005	0.001	1
25	8	S	221	0.002	1.122 *	1
26	9	S	232	-0.002	0.000	1
27	10	S	210	0.002	0.001	1
28	11	K	210 spiked	0.146	0.142	1
29	908	QC1	CRM N03 [0.030] Old	0.010 F	-0.001	1
30	909	QC2	CRM N03 [0.030] New	0.030	0.001	1
31	0	B	Baseline	-0.002	-0.002	1
32	12	S	331	0.004	0.001	1
33	13	S	342	0.005	0.001	1
34	14	S	397	0.001	0.001	1
35	15	S	320	0.003	0.001	1
36	16	S	364	0.001	0.000	1
37	17	S	386	0.003	0.000	1
38	18	S	463	0.006	0.001	1
39	19	S	353	0.000	0.000	1
40	20	S	452	0.005	0.002	1
41	21	S	408	-0.002	0.001	1
42	22	K	408 spiked	0.013	0.139	1

43	908	QC1	CRM N03 [0.030] Old	0.011 F	-0.001	1
44	909	QC2	CRM N03 [0.030] New	0.030	0.001	1
45	0	B	Baseline	-0.002	-0.002	1
46	23	S	419	-0.002	0.001	1
47	24	S	430	-0.003	0.000	1
48	25	S	441	0.000	0.001	1
49	26	S	474	0.001	0.052	1
50	27	S	518	0.001	0.001	1
51	28	S	529	0.001	0.001	1
52	29	S	540	0.001	0.003	1
53	30	S	507	-0.001	0.001	1
54	31	S	562	0.002	0.000	1
55	32	S	551	0.002	0.001	1
56	33	K	551 spiked	0.129	0.126	1
57	908	QC1	CRM N03 [0.030] Old	0.010 F	-0.001	1
58	909	QC2	CRM N03 [0.030] New	0.030	0.001	1
59	0	B	Baseline	-0.002	-0.002	1
60	34	S	573	-0.003	0.000	1
61	35	S	584	0.002	0.001	1
62	36	S	485	0.002	0.001	1
63	37	S	617	0.001	0.001	1
64	38	S	375	0.003	0.002	1
65	39	S	595	0.001	0.001	1
66	40	K	595 spiked	0.123	0.135	1
67	908	QC1	CRM N03 [0.030] Old	0.010 F	-0.001	1
68	909	QC2	CRM N03 [0.030] New	0.029	0.002	1
69	0	B	Baseline	-0.002	-0.002	1
70	921	D	Drift	0.199	0.200	1
71	0	B	Final Base	-0.002	-0.002	1

Corrections

Channel	1	3
Baseline	done	done
Drift	none	none
Carryover	done	done
%	0.54	0.17

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.138	0.137
Recovery %	92.267	91.509
210 spiked Added	0.150	0.150
210 spiked Found	0.145	0.141
Recovery %	96.446	94.002
408 spiked Added	0.150	0.150
408 spiked Found	0.015	0.138
Recovery %	9.942	92.286
551 spiked Added	0.150	0.150
551 spiked Found	0.126	0.125
Recovery %	84.215	83.554
595 spiked Added	0.150	0.150
595 spiked Found	0.122	0.134
Recovery %	81.432	89.567

Calibration

Channel 1 NOx

Type	Observed	Calculated	Target	% Error
1C	57943	0.200	0.200	-0.24
2C	31145	0.101	0.100	0.70
3C	17521	0.050	0.050	0.91
4C	10679	0.025	0.025	0.90
5C	6603	0.010	0.010	1.94
6C	5031	0.004	0.005	-12.06
7C	3708	0.000	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 1.0000
 Carryover (%): 0.5440
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 3.838643E+003
 b = 2.711782E+005

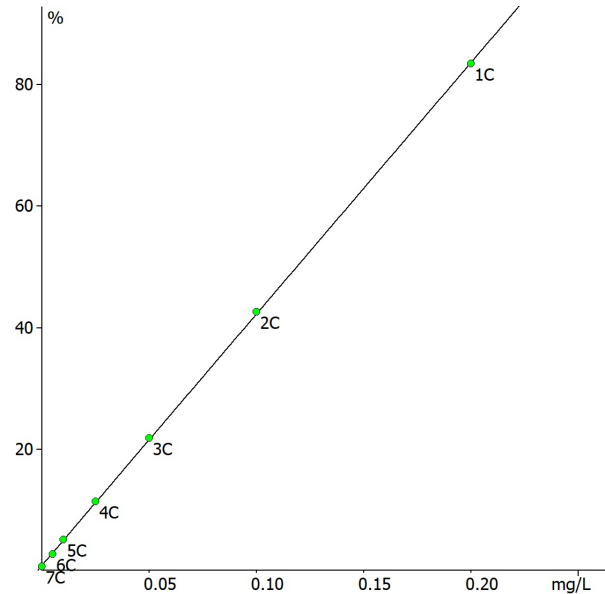

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	57945	0.200	0.200	-0.14
2C	31184	0.101	0.100	0.99
3C	17184	0.049	0.050	-1.31
4C	10394	0.024	0.025	-2.82
5C	6690	0.011	0.010	6.29
6C	5167	0.005	0.005	0.21
7C	3810	0.000	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 1.0000
 Carryover (%): 0.1690
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 3.808898E+003
 b = 2.710557E+005

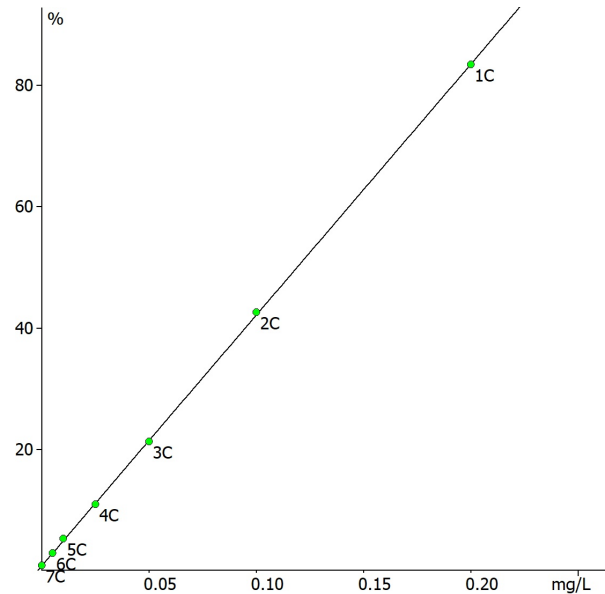

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero