

General

Name of Run	240305_NP1_QCSettings_3CNull.run	Name of Analysis	NOx_P04.ANL
Date of Report	3/6/2024	System	AA500
Run Start	3/5/2024 6:44:36 PM	Run Stop	9:05:27 PM
Operator Run/Rec.	CLMR / CLMR	Software Version	AACE 8.02 SP2
Comment	Recalculate from Run240305_NP1_QCSettings.run te from Run240305_NP1.run		

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9999	0.9999	
AUFS/Gain	0.40/25	0.10/105	Data
Offset	320	5005	Field
Sensitivity	0.3221	0.0796	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.002	0.005	0
1	921	P	Primer	0.200	0.198	1
2	921	D	Drift	0.199	0.199	1
3	921	C	0.200	0.200	0.199	1
4	901	C	0.100	0.100	0.101	1
5	0	XC	Cal. (3)	0.096 N	0.051	1
6	903	C	0.025	0.025	0.027	1
7	904	C	0.010	0.011	0.009	1
8	905	C	0.005	0.006	0.004	1
9	906	C	Blank	-0.002	-0.001	1
10	0	B	Baseline	-0.002	0.005	1
11	921	H1	High	0.200	0.198	1
12	905	L1	Low	0.006	0.004	1
13	905	L1	Low	0.006	0.004	1
14	906	S	Blank	-0.001	-0.001	1
15	907	K	Blank spiked	0.144	0.144	1
16	908	QC1	CRM N03	0.133	0.029	1
17	909	QC2	CRM P04	0.036	0.008	1
18	1	S	31	-0.003	0.005	1
19	2	S	42	0.007	0.001	1
20	3	S	75	0.014	0.001	1
21	4	S	282	0.004	0.003	1
22	5	S	64	-0.002	0.001	1
23	6	S	20	0.004	0.001	1
24	7	S	9	0.005	0.001	1
25	8	S	53	-0.001	0.001	1
26	9	S	152	0.002	0.001	1
27	10	S	97	-0.002	0.001	1
28	11	K	97 spiked	0.044	0.150	1
29	908	QC1	CRM N03	0.133	0.027	1
30	909	QC2	CRM P04	0.037	0.007	1
31	0	B	Baseline	-0.002	0.005	1
32	12	S	126	0.004	0.002	1
33	13	S	334	0.010	0.001	1
34	14	S	243	0.005	0.000	1
35	15	S	139	-0.002	0.001	1
36	16	S	217	0.006	0.001	1
37	17	S	191	0.007	0.001	1
38	18	S	256	0.006	0.001	1
39	19	S	204	-0.003	0.000	1
40	20	S	269	0.003	0.000	1
41	21	K	269 spiked	0.139	0.147	1
42	908	QC1	CRM N03	0.133	0.027	1

43	909	QC2	CRM P04	0.037	0.008	1
44	0	B	Baseline	-0.002	0.005	1
45	22	S	321	0.001	0.001	1
46	23	S	308	0.010	0.001	1
47	24	S	230	0.008	0.001	1
48	25	S	360	0.009	0.001	1
49	26	S	165	0.008	0.000	1
50	27	S	386	0.009	0.001	1
51	28	S	399	0.004	0.003	1
52	29	S	425	0.126	0.003	1
53	30	S	373	0.004	0.000	1
54	31	S	347	0.034	0.001	1
55	32	K	347 spiked	0.178	0.153	1
56	908	QC1	CRM N03	0.133	0.028	1
57	909	QC2	CRM P04	0.037	0.007	1
58	0	B	Baseline	-0.002	0.005	1
59	33	S	438	0.004	0.005	1
60	34	S	178	0.001	0.001	1
61	35	S	516	0.008	0.005	1
62	36	S	503	0.006	0.003	1
63	37	S	451	0.006	0.003	1
64	38	S	477	0.024	0.000	1
65	39	S	555	0.049	0.002	1
66	40	S	412	0.006	0.005	1
67	41	S	529	0.004	0.004	1
68	42	S	568	0.010	0.000	1
69	43	K	568 spiked	0.168	0.152	1
70	908	QC1	CRM N03	0.134	0.027	1
71	909	QC2	CRM P04	0.036	0.007	1
72	0	B	Baseline	-0.002	0.005	1
73	44	S	633	0.004	0.001	1
74	45	S	594	0.006	0.001	1
75	46	S	620	0.007	0.002	1
76	47	S	542	0.008	0.002	1
77	48	S	295	0.003	0.002	1
78	49	S	490	0.011	0.003	1
79	50	S	581	-0.004	0.003	1
80	51	S	86	0.005	0.005	1
81	52	K	86 spiked	0.140	0.157	1
82	908	QC1	CRM N03	0.135	0.032	1
83	909	QC2	CRM P04	0.035	0.012	1
84	0	B	Baseline	-0.002	0.005	1
85	910	S	Blank	0.002	0.005	1
86	910	S	Blank	0.002	0.005	1
87	910	S	Blank	0.001	0.005	1
88	910	S	Blank	0.001	0.004	1
89	910	S	Blank	0.001	0.004	1
90	910	S	Blank	0.000	0.004	1
91	910	S	Blank	0.001	0.004	1
92	910	S	Blank	0.000	0.004	1
93	910	S	Blank	0.001	0.004	1
94	910	S	Blank	0.000	0.004	1
95	908	QC1	CRM N03	0.135	0.034	1
96	909	QC2	CRM P04	0.035	0.013	1
97	0	B	Baseline	-0.002	0.005	1
98	921	D	Drift	0.202	0.205	1
99	0	B	Final Base	-0.002	0.005	1

Corrections

Channel	1	3
Baseline	done	done
Drift	none	none
Carryover	done	done
%	2.98	0.34

Spike Recovery

Blank spiked Added	0.150	0.150
Blank spiked Found	0.145	0.145
Recovery %	96.838	96.681
97 spiked Added	0.150	0.150
97 spiked Found	0.046	0.149
Recovery %	30.509	99.284
269 spiked Added	0.150	0.150


269 spiked Found	0.136	0.147
Recovery %	90.681	97.865
347 spiked Added	0.149	0.150
347 spiked Found	0.145	0.152
Recovery %	97.131	101.477
568 spiked Added	0.150	0.150
568 spiked Found	0.158	0.152
Recovery %	105.452	101.196
86 spiked Added	0.150	0.150
86 spiked Found	0.135	0.152
Recovery %	90.324	101.608

Calibration

Channel 1 NOx

Type	Observed	Calculated	Target	% Error
1C	56103	0.200	0.200	0.06
2C	29802	0.100	0.100	-0.43
4C	10397	0.025	0.025	1.52
5C	6519	0.011	0.010	5.53
6C	5332	0.006	0.005	20.31
7C	3332	-0.002	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 0.9999
 Carryover (%): 2.9829
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = $3.758540E+003$
 b = $2.615699E+005$

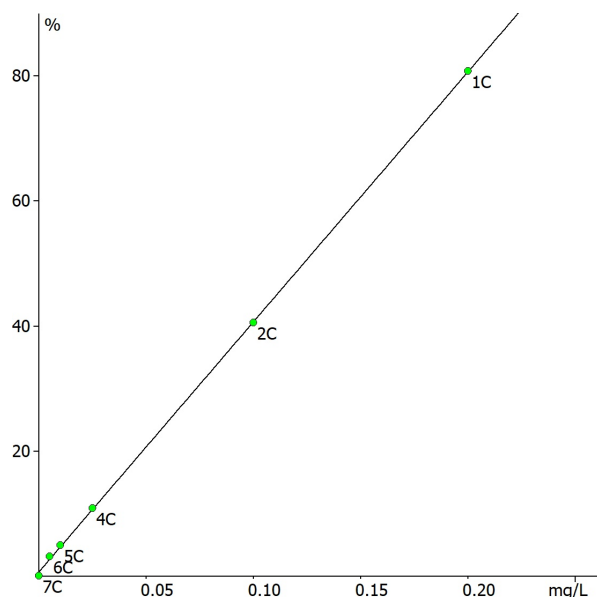

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58315	0.199	0.200	-0.42
2C	30431	0.101	0.100	0.81
3C	16339	0.051	0.050	2.22
4C	9396	0.027	0.025	6.47
5C	4519	0.009	0.010	-5.85
6C	2983	0.004	0.005	-20.05
7C	1532	-0.001	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9999
 Carryover (%): 0.3367
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = $1.849682E+003$
 b = $2.835057E+005$

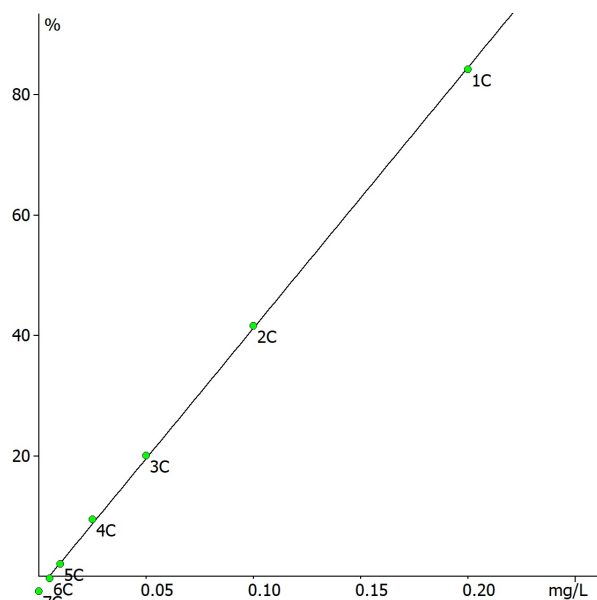

Legend

* ... Sample offscale
+ ... Result higher than sample limit
- ... Result lower than sample limit
P ... Standard passed
F ... Standard failed
N ... Value not calculated or not used
R ... Resample after offscale
M ... Peak marker moved manually
d ... Pre-Diluted sample
D ... Post-Diluted sample
H ... Dual calibration used
Z ... Report negative values as zero