

General

Name of Run	240307_Q3_TKN1R1.RUN	Name of Analysis	TKN.ANL
Date of Report	3/7/2024	System	AA500
Run Start	3/7/2024 5:07:41 PM	Run Stop	7:41:29 PM
Operator Run/Rec. CLMR	/	Software Version	AACE 8.02 SP2
Comment			

Parameters

Channel	1		
Method	TKN	Dilution	
Unit	mg/L		
Calib. Fit	Linear		
Corr.Coeff.(r)	0.9958		
AUFS/Gain	0.18/56	Data	
Offset	-1982	Field	
Sensitivity	0.1136	'#2	
Sample Limit 1	-		
Sample Limit 2	-		

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.138	0	
1	921	P	Primer	2.133	1	
2	921	D	Drift	2.174	1	
3	921	C	2.22	2.144	1	
4	901	C	1.66	1.713	1	
5	902	C	1.11	1.165	1	
6	903	C	0.55	0.634	1	
7	904	C	0.22	0.105	1	
8	905	C	Blank	-0.001	1	
9	0	B	Baseline	-0.138	1	
10	921	H1	High	2.104	1	
11	904	L1	Low	0.090	1	
12	904	L1	Low	0.098	1	
13	905	S	Blank	-0.027	1	
14	906	K	Blank spiked (No Acid St~	0.918	1	
15	905	S	Blank	-0.068	1	
16	907	K	Blank spiked (Acid Stock~	0.973	1	
17	908	QC1	CRM [1.5 mg/L]	1.350	1	
18	1	S	426	0.124	1	
19	2	S	231	0.038	1	
20	3	S	361	-0.017	1	
21	4	S	374	-0.065	1	
22	5	S	387	0.028	1	
23	6	S	270	-0.035	1	
24	7	S	296	0.510	1	
25	8	S	309	0.230	1	
26	9	S	400	-0.100	1	
27	10	S	322	0.412	1	
28	11	K	322 spiked	1.240	1	
29	908	QC1	CRM [1.5 mg/L]	1.443	1	
30	0	B	Baseline	-0.138	1	
31	12	S	166	0.007	1	
32	13	S	65	-0.019	1	
33	14	S	192	-0.088	1	
34	15	S	76	0.106	1	
35	16	S	205	-0.152	1	
36	17	S	43	-0.011	1	
37	18	S	257	-0.061	1	
38	19	S	32	-0.009	1	
39	20	S	127	-0.138	1	
40	21	S	283	-0.070	1	
41	22	K	283 spiked	0.893	1	
42	908	QC2	CRM [1.5 mg/L]	1.460	1	

43	0	B	Baseline	-0.138 M	1
44	23	S	153	0.021	1
45	24	S	87	-0.050	1
46	25	S	54	-0.078	1
47	26	S	140	-0.103	1
48	27	S	10	0.083	1
49	28	S	244	0.021	1
50	29	S	218	-0.032	1
51	30	S	335	0.729	1
52	31	S	21	0.082	1
53	32	S	98	-0.009	1
54	33	K	98 spiked	0.928	1
55	908	QC1	CRM [1.5 mg/L]	1.426	1
56	0	B	Baseline	-0.138 M	1
57	34	S	179	0.456	1
58	35	S	348	0.160	1
59	36	S	439	0.308	1
60	37	S	491	0.304	1
61	38	S	569	0.355	1
62	39	S	517	0.148	1
63	40	S	452	0.044	1
64	41	S	504	0.109	1
65	42	S	478	-0.003	1
66	43	S	556	0.116	1
67	44	K	556 spiked	1.137	1
68	908	QC1	CRM [1.5 mg/L]	1.290	1
69	0	B	Baseline	-0.138	1
70	45	S	621	0.226	1
71	46	S	634	0.057	1
72	47	S	582	0.197	1
73	48	S	530	-0.031	1
74	49	S	413	0.078	1
75	50	S	595	0.264	1
76	51	S	543	-0.011	1
77	52	K	543 spiked	0.984	1
78	908	QC1	CRM [1.5 mg/L]	1.244	1
79	0	B	Baseline	-0.138	1
80	922	S	Blank	-0.081	1
81	922	S	Blank	-0.114	1
82	922	S	Blank	-0.137	1
83	922	S	Blank	-0.179	1
84	922	S	Blank	-0.135	1
85	922	S	Blank	-0.128	1
86	922	S	Blank	-0.164	1
87	922	S	Blank	-0.124	1
88	922	S	Blank	-0.146	1
89	922	S	Blank	-0.189	1
90	908	QC1	CRM [1.5 mg/L]	1.282	1
91	0	B	Baseline	-0.138	1
92	921	D	Drift	2.174	1
93	0	B	Final Base	-0.138	1

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	none
%	-0.37

Spike Recovery

Blank spiked (No Acid Stock B) Added	1.001
Blank spiked (No Acid Stock B) Found	0.945
Recovery %	94.367
Blank spiked (Acid Stock B) Added	1.010
Blank spiked (Acid Stock B) Found	1.041
Recovery %	103.098
322 spiked Added	0.998
322 spiked Found	0.829
Recovery %	83.034
283 spiked Added	1.001
283 spiked Found	0.963
Recovery %	96.155


98 spiked Added	1.001
98 spiked Found	0.937
Recovery %	93.623
556 spiked Added	1.000
556 spiked Found	1.021
Recovery %	102.106
543 spiked Added	1.001
543 spiked Found	0.995
Recovery %	99.444

Calibration

Channel 1	TKN
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	45584	2.144	2.220	-3.44
2C	37671	1.713	1.660	3.17
3C	27617	1.165	1.110	4.96
4C	17875	0.634	0.550	15.35
5C	8161	0.105	0.220	-52.13
6C	6210	-0.001	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9958
 Carryover (%): no
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 6.227417E+003
 b = 1.835998E+004

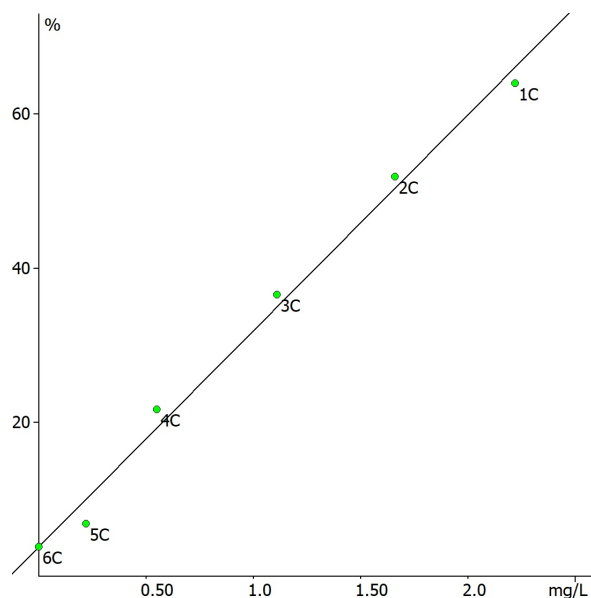

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero