

General

Name of Run	240430_Q4_NP1R2.RUN	Name of Analysis	NOx_P04.ANL
Date of Report	5/1/2024	System	AA500
Run Start	4/30/2024 6:09:56 PM	Run Stop	7:44:22 PM
Operator Run/Rec. CLMR	/	Software Version	AACE 8.02 SP2
Comment			

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9998	0.9999	
AUFS/Gain	0.42/24	0.09/107	Data
Offset	203	4816	Field
Sensitivity	0.3331	0.0782	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.004	0.000	0
1	921	P	Primer	0.198	0.198	1
2	921	D	Drift	0.203	0.199	1
3	921	C	0.200	0.200	0.200	1
4	901	C	0.100	0.099	0.101	1
5	902	C	0.050	0.050	0.049	1
6	903	C	0.025	0.027	0.024	1
7	904	C	0.010	0.009	0.009	1
8	905	C	0.005	0.007	0.005	1
9	906	C	Blank	-0.002	0.002 M	1
10	0	B	Baseline	-0.004	0.000	1
11	921	H1	High	0.200	0.198	1
12	905	L1	Low	0.007	0.005	1
13	905	L1	Low	0.007	0.005	1
14	906	S	Blank	-0.002	0.033	1
15	907	K	Blank spiked	0.144	0.148	1
16	908	QC1	CRM [0.15/0.03]	0.148 F	0.039	1
17	909	QC1	CRM Old	0.153 F	0.119	1
18	1	S	35	0.006	0.013	1
19	2	S	191	0.010	0.008	1
20	3	S	100	0.009	0.007	1
21	4	S	22	0.009	0.014	1
22	5	S	74	0.021	0.039	1
23	6	S	48	0.010	0.009	1
24	7	S	9	0.010	0.000	1
25	8	S	126	0.015	0.004	1
26	9	S	139	0.015	0.005	1
27	10	S	152	0.013	0.003	1
28	11	K	152 spiked	0.158	0.149	1
29	908	QC1	CRM [0.15/0.03]	0.146 F	0.032	1
30	909	QC1	CRM Old	0.152 F	0.121	1
31	0	B	Baseline	-0.004	0.000	1
32	12	S	165	0.004	0.001	1
33	13	S	230	0.005	0.006	1
34	14	S	373	0.003	0.017	1
35	15	S	243	0.013	0.025	1
36	16	S	217	0.013	0.005	1
37	17	S	269	0.009	0.003	1
38	18	S	282	-0.003	0.000	1
39	19	S	256	0.001	0.006	1
40	20	S	321	0.002	0.009	1
41	21	S	87	0.005	0.005	1
42	22	K	87 spiked	0.148	0.151	1

43	908	QC1	CRM [0.15/0.03]	0.148 F	0.032	1
44	909	QC1	CRM Old	0.152 F	0.120	1
45	0	B	Baseline	-0.004	0.000	1
46	23	S	113	0.007	0.001	1
47	24	S	308	0.001	0.008	1
48	25	S	48 Dup	0.013	0.024	1
49	26	S	61	0.004	0.004	1
50	27	K	61 spiked	0.149	0.149	1
51	908	QC1	CRM [0.15/0.03]	0.148 F	0.030	1
52	909	QC1	CRM Old	0.152 F	0.122	1
53	922	S	Blank	0.000	0.006	1
54	922	S	Blank	0.000	0.003	1
55	922	S	Blank	0.000	0.001	1
56	922	S	Blank	0.000	0.002	1
57	922	S	Blank	0.000	0.002	1
58	922	S	Blank	0.000	0.002	1
59	922	S	Blank	0.000	0.002	1
60	922	S	Blank	0.000	0.004	1
61	922	S	Blank	-0.001 M	0.004	1
62	922	S	Blank	0.000	0.004	1
63	0	B	Baseline	-0.004	0.000	1
64	921	D	Drift	0.198	0.199	1
65	0	B	Final Base	-0.004	0.000	1

Corrections

Channel	1	3
Baseline	done	done
Drift	none	none
Carryover	done	done
%	0.47	0.08

Spike Recovery


Blank spiked Added	0.150	0.149
Blank spiked Found	0.146	0.115
Recovery %	96.976	76.988
152 spiked Added	0.150	0.150
152 spiked Found	0.145	0.147
Recovery %	96.690	97.842
87 spiked Added	0.150	0.150
87 spiked Found	0.143	0.146
Recovery %	95.577	97.637
61 spiked Added	0.150	0.150
61 spiked Found	0.145	0.145
Recovery %	96.545	96.855

Calibration

Channel 1	NOx
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	56261	0.200	0.200	0.02
2C	30097	0.099	0.100	-0.56
3C	17307	0.050	0.050	0.54
4C	11147	0.027	0.025	6.35
5C	6463	0.009	0.010	-14.22
6C	6053	0.007	0.005	40.04
7C	3734	-0.002	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 0.9998
 Carryover (%): 0.4718
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 4.231732E+003
 b = 2.601022E+005

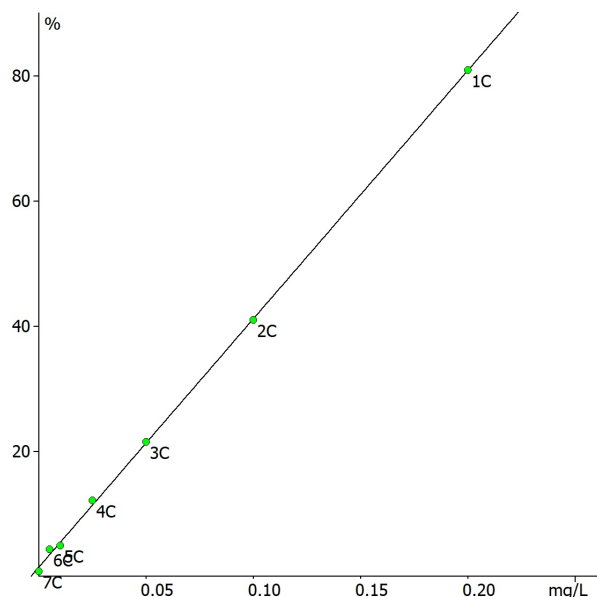

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58604	0.200	0.200	0.01
2C	31050	0.101	0.100	0.69
3C	16790	0.049	0.050	-1.44
4C	9765	0.024	0.025	-4.17
5C	5596	0.009	0.010	-10.70
6C	4609	0.005	0.005	7.45
7C	3607	0.002	0.000	-

Calibration Fit: Linear
Correlation Coefficient: 0.9999
Carryover (%): 0.0800
Calibration Equation: $y = bx + a$
y = Peak Height Digit. Units
x = Concentration Unit mg/L
a = 3.118652E+003
b = 2.774114E+005

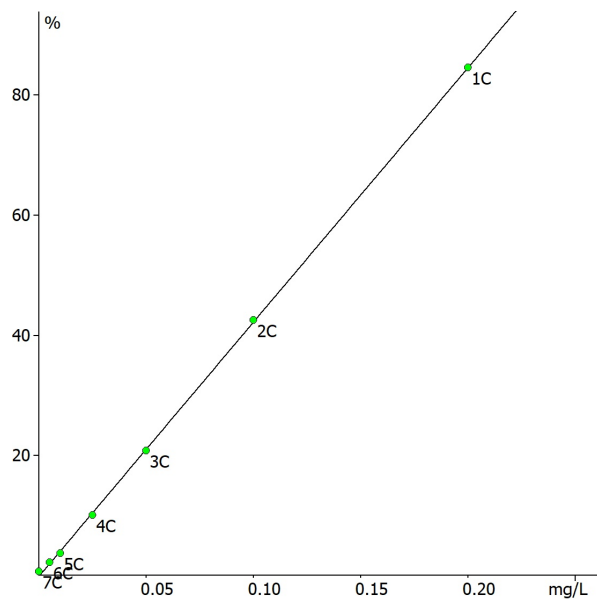

Legend

* ... Sample offscale
+ ... Result higher than sample limit
- ... Result lower than sample limit
P ... Standard passed
F ... Standard failed
N ... Value not calculated or not used
R ... Resample after offscale
M ... Peak marker moved manually
d ... Pre-Diluted sample
D ... Post-Diluted sample
H ... Dual calibration used
Z ... Report negative values as zero