

General

Name of Run	240506_Q4_TKN1.run	Name of Analysis	TKN.ANL
Date of Report	5/10/2024	System	AA500
Run Start	5/6/2024 4:15:43 PM	Run Stop	6:18:30 PM
Operator Run/Rec. CLMR		Software Version	AACE 8.02 SP2
Comment			

Parameters

Channel	1		
Method	TKN	Dilution	
Unit	mg/L		
Calib. Fit	Linear		
Corr.Coeff.(r)	0.9828		
AUFS/Gain	0.16/61	Data	
Offset	-1279	Field	
Sensitivity	0.1381	'#2	
Sample Limit 1	-		
Sample Limit 2	-		

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.240	0	
1	921	P	Primer	2.367	1	
2	921	D	Drift	2.365	1	
3	921	C	2.22	2.373	1	
4	901	C	1.66	1.508	1	
5	902	C	1.11	0.923	1	
6	903	C	0.55	0.726	1	
7	904	C	0.22	0.324	1	
8	905	C	Blank	-0.094	1	
9	0	B	Baseline	-0.240	1	
10	921	H1	High	2.380	1	
11	904	L1	Low	0.320	1	
12	904	L1	Low	0.320	1	
13	905	S	Blank	-0.089	1	
14	906	K	Blank spiked (No Acid St~	0.906	1	
15	907	QC1	CRM [1.5 mg/L]	1.477	1	
16	1	S	218	-0.008	1	
17	2	S	452	0.076	1	
18	3	S	257	-0.024	1	
19	4	S	374	0.080	1	
20	5	S	166	-0.036	1	
21	6	S	10	0.008	1	
22	7	S	101	0.109	1	
23	8	S	322	0.063	1	
24	9	S	62	0.961	1	
25	10	S	36 SpTest NoAcid	1.041	1	
26	11	S	36 SpTest Acid	1.037	1	
27	12	S	36	0.051	1	
28	13	K	36 spiked	0.864	1	
29	907	QC1	CRM [1.5 mg/L]	1.469	1	
30	0	B	Baseline	-0.240	1	
31	14	S	127	-0.036	1	
32	15	S	140	0.769	1	
33	16	S	335	0.343	1	
34	17	S	75	0.144	1	
35	18	S	244	0.128	1	
36	19	S	205	0.252	1	
37	20	S	296	0.134	1	
38	21	S	49	0.329	1	
39	22	S	192	0.192	1	
40	23	S	114 SpTest NoAcid	1.064	1	
41	24	S	114 SpTest Acid	1.032	1	
42	25	S	114	0.040	1	

43	26	K	114 spiked	1.023	1
44	907	QC2	CRM [1.5 mg/L]	1.473	1
45	0	B	Baseline	-0.240	1
46	27	S	361	0.165	1
47	28	S	153	0.303	1
48	29	S	179	0.040	1
49	30	S	23	0.002	1
50	31	S	478	0.129	1
51	32	S	88	0.278	1
52	33	S	309	-0.003	1
53	34	S	270	-0.051	1
54	35	S	348	0.049	1
55	36	S	465	0.044	1
56	37	S	283	0.218	1
57	38	S	231 SpTest NoAcid	1.033	1
58	39	S	231 SpTest Acid	0.997	1
59	40	S	231	0.016	1
60	41	K	231 spiked	0.835	1
61	907	QC1	CRM [1.5 mg/L]	1.466	1
62	0	B	Baseline	-0.240	1
63	905	S	Blank	-0.096	1
64	905	S	Blank	-0.096	1
65	905	S	Blank	-0.097	1
66	905	S	Blank	-0.094	1
67	905	S	Blank	-0.096	1
68	905	S	Blank	-0.096	1
69	905	S	Blank	-0.090	1
70	905	S	Blank	-0.091	1
71	905	S	Blank	-0.097	1
72	905	S	Blank	-0.098	1
73	907	QC1	CRM [1.5 mg/L]	1.468	1
74	0	B	Baseline	-0.240	1
75	921	D	Drift	2.365	1
76	0	B	Final Base	-0.240	1

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	done
%	0.33

Spike Recovery


Blank spiked (No Acid Stock B) Added	1.001
Blank spiked (No Acid Stock B) Found	0.996
Recovery %	99.431
36 spiked Added	1.000
36 spiked Found	0.813
Recovery %	81.262
114 spiked Added	1.001
114 spiked Found	0.983
Recovery %	98.211
231 spiked Added	1.001
231 spiked Found	0.819
Recovery %	81.836

Calibration

Channel 1 TKN

Type	Observed	Calculated	Target	% Error
1C	58664	2.373	2.220	6.91
2C	40337	1.508	1.660	-9.16
3C	27942	0.923	1.110	-16.87
4C	23777	0.726	0.550	32.02
5C	15256	0.324	0.220	47.18
6C	6407	-0.094	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9828
 Carryover (%): 0.3291
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 8.398477E+003
 b = 2.117919E+004

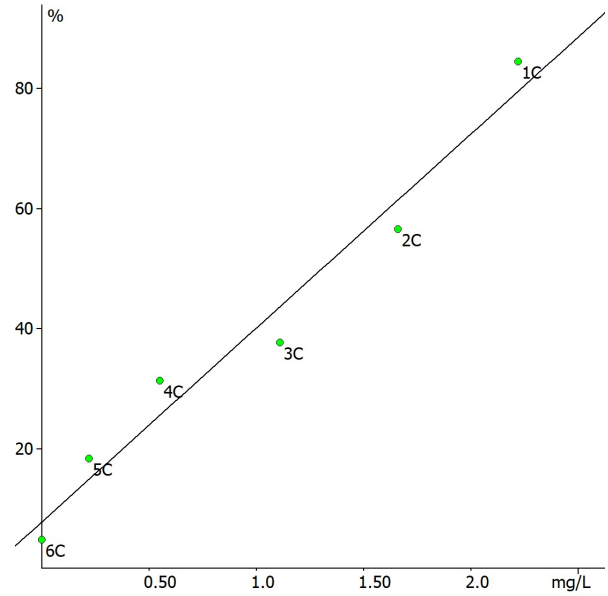

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero