

General

Name of Run	240509_Q4_TKN2.run	Name of Analysis	TKN.ANL
Date of Report	5/10/2024	System	AA500
Run Start	5/9/2024 1:57:07 PM	Run Stop	4:48:34 PM
Operator Run/Rec. CLMR		Software Version	AACE 8.02 SP2
Comment			

Parameters

Channel	1		
Method	TKN	Dilution	
Unit	mg/L		
Calib. Fit	Linear		
Corr.Coeff.(r)	0.9932		
AUFS/Gain	0.14/69	Data	
Offset	-1797	Field	
Sensitivity	0.1104	'#2	
Sample Limit 1	-		
Sample Limit 2	-		

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.188	0	
1	921	P	Primer	2.107	1	
2	921	D	Drift	2.082	1	
3	921	C	2.22	2.088	1	
4	901	C	1.66	1.841	1	
5	902	C	1.11	1.110	1	
6	903	C	0.55	0.537	1	
7	904	C	0.22	0.218	1	
8	905	C	Blank	-0.033	1	
9	0	B	Baseline	-0.188	1	
10	921	H1	High	2.083	1	
11	904	L1	Low	0.243	1	
12	904	L1	Low	0.243	1	
13	905	S	Blank	-0.028	1	
14	906	K	Blank Spike Dig.	1.009	1	
15	907	QC1	CRM [1.5 mg/L]	1.638	1	
16	1	S	387	0.098	1	
17	2	S	530	0.107	1	
18	3	S	608	0.041	1	
19	4	S	426	0.128	1	
20	5	S	517	-0.064	1	
21	6	K	517 spiked	1.117	1	
22	7	S	621	-0.012	1	
23	8	S	634	0.067	1	
24	9	S	582	0.367	1	
25	10	S	595	-0.104	1	
26	11	S	413	0.024	1	
27	12	K	413 spiked	1.148	1	
28	907	QC1	CRM [1.5 mg/L]	1.619	1	
29	0	B	Baseline	-0.188	1	
30	13	S	439	-0.060	1	
31	14	S	543	0.422	1	
32	15	S	504	0.064	1	
33	16	S	400	0.003	1	
34	17	S	569	0.051	1	
35	18	K	569 spiked	1.174	1	
36	907	QC1	CRM [1.5 mg/L]	1.631	1	
37	19	S	C. 2.22	2.680 *	1	
38	20	S	C. 1.66	1.658	1	
39	21	S	C. 1.11	1.052	1	
40	22	S	C. 0.55	0.845	1	
41	23	S	C. 0.22	0.397	1	
42	24	S	Blnk. Spike	1.026	1	

43	907	QC1	CRM [1.5 mg/L]	1.633	1
44	0	B	Baseline	-0.188	1
45	25	S	218	0.076	1
46	26	S	452	0.196	1
47	27	S	257	0.076	1
48	28	S	374	0.197	1
49	29	S	166	-0.062	1
50	30	S	10	0.041	1
51	31	S	101	0.137	1
52	32	S	322	0.076	1
53	33	S	62	1.061	1
54	34	S	36 SpTest NoAcid	1.138	1
55	35	S	36 SpTest Acid	1.146	1
56	36	S	36	0.079	1
57	37	K	36 spiked	0.989	1
58	907	QC1	CRM [1.5 mg/L]	1.618	1
59	0	B	Baseline	-0.188	1
60	38	S	127	0.025	1
61	39	S	140	0.906	1
62	40	S	335	0.464	1
63	41	S	75	0.228	1
64	42	S	244	0.231	1
65	43	S	205	0.375	1
66	44	S	296	0.251	1
67	45	S	49	0.452	1
68	46	S	192	0.317	1
69	47	S	114 SpTest NoAcid	1.249	1
70	48	S	114 SpTest Acid	1.181	1
71	49	S	114	0.108	1
72	50	K	114 spiked	1.179	1
73	907	QC2	CRM [1.5 mg/L]	1.656	1
74	0	B	Baseline	-0.188	1
75	51	S	361	0.283	1
76	52	S	153	0.441	1
77	53	S	179	-0.080	1
78	54	S	23	0.106	1
79	55	S	478	0.229	1
80	56	S	88	0.423	1
81	57	S	309	0.094	1
82	58	S	270	0.061	1
83	59	S	348	0.156	1
84	60	S	465	0.160	1
85	911	S	283	0.337	1
86	912	S	231 SpTest NoAcid	1.209	1
87	913	S	231 SpTest Acid	1.131	1
88	914	S	231	0.925	1
89	915	K	231 spiked	-0.390	1
90	907	QC1	CRM [1.5 mg/L]	1.639	1
91	0	B	Baseline	-0.188	1
92	905	S	Blank	-0.019	1
93	905	S	Blank	-0.034	1
94	905	S	Blank	0.007	1
95	905	S	Blank	0.002	1
96	905	S	Blank	0.048	1
97	905	S	Blank	0.027	1
98	905	S	Blank	0.017	1
99	905	S	Blank	-0.012	1
100	905	S	Blank	-0.008	1
101	905	S	Blank	-0.390	1
102	907	QC1	CRM [1.5 mg/L]	1.618	1
103	0	B	Baseline	-0.188	1
104	921	D	Drift	2.082	1
105	0	B	Final Base	-0.188	1

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	done
%	1.35

Spike Recovery

Blank Spike Dig. Added	1.001
------------------------	-------


Blank Spike Dig. Found	1.037
Recovery %	103.600
517 spiked Added	1.001
517 spiked Found	1.180
Recovery %	117.883
413 spiked Added	1.001
413 spiked Found	1.124
Recovery %	112.369
569 spiked Added	1.000
569 spiked Found	1.123
Recovery %	112.249
36 spiked Added	1.000
36 spiked Found	0.910
Recovery %	90.998
114 spiked Added	1.000
114 spiked Found	1.072
Recovery %	107.187
231 spiked Added	0.994
231 spiked Found	-1.315
Recovery %	-132.277

Calibration

Channel 1 TKN

Type	Observed	Calculated	Target	% Error
1C	53989	2.088	2.220	-5.93
2C	48611	1.841	1.660	10.91
3C	32699	1.110	1.110	-0.04
4C	20237	0.537	0.550	-2.45
5C	13303	0.218	0.220	-1.04
6C	7845	-0.033	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9932
 Carryover (%): 1.3467
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 8.567937E+003
 b = 2.174937E+004

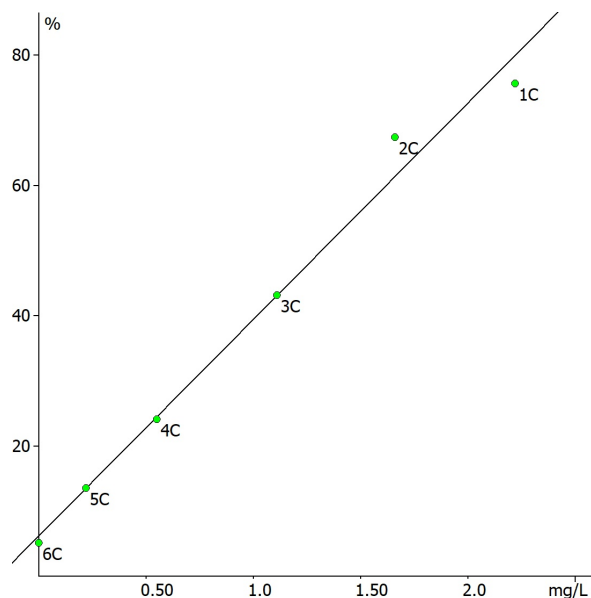

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero