

General

Name of Run	240514_Q4_NP2R1_nobblankP.run	Name of Analysis	NOx_P04.ANL
Date of Report	5/14/2024	System	AA500
Run Start	5/14/2024 3:40:26 PM	Run Stop	5:14:56 PM
Operator Run/Rec.	CLMR / tac	Software Version	AACE 8.02 SP2
Comment	Recalculate from Run240514_Q4_NP2.run		

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9994	0.9999	
AUFS/Gain	0.42/24	0.09/107	Data
Offset	194	4832	Field
Sensitivity	0.3294	0.0774	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	0.001	0.001	0
1	922	P	Primer	0.200	0.198	1
2	922	D	Drift	0.201	0.199	1
3	922	C	0.200	0.201	0.199	1
4	911	C	0.100	0.098	0.102	1
5	912	C	0.050	0.050	0.050	1
6	913	C	0.025	0.022	0.024	1
7	914	C	0.010	0.011	0.010	1
8	915	C	0.005	0.004	0.005	1
9	0	XC	Cal. (1)	0.004	0.030 N	1
10	0	B	Baseline	0.001	0.001	1
11	922	H1	High	0.201	0.197	1
12	915	L1	Low	0.004	0.003	1
13	915	L1	Low	0.004	0.003	1
14	916	S	Blank	0.004	0.033	1
15	917	K	Blank spiked	0.151	0.154	1
16	918	QC1	CRM [0.15/0.03]	0.150	0.031	1
17	1	S	178	0.019	-0.001	1
18	2	S	386	0.014	-0.002	1
19	3	S	347	0.012	-0.001	1
20	4	S	204	0.014	-0.001	1
21	5	S	568	0.016	-0.001	1
22	6	S	633	0.014	-0.001	1
23	7	S	503	0.021	-0.001	1
24	8	S	464	0.012	-0.001	1
25	9	S	607	0.013	-0.001	1
26	10	S	425	0.015	0.001	1
27	11	K	425 spiked	0.161	0.156	1
28	918	QC1	CRM [0.15/0.03]	0.150	0.038	1
29	0	B	Baseline	0.001	0.001	1
30	12	S	516	0.014	0.000	1
31	13	S	477	0.012	0.000	1
32	14	S	295	0.017	0.000	1
33	15	S	581	0.016	0.000	1
34	16	S	438	0.018	0.000	1
35	17	S	412	0.014	0.002	1
36	18	S	594	0.016	0.000	1
37	19	S	529	0.012	0.000	1
38	20	S	620	0.016	0.001	1
39	21	S	542	0.054	0.020	1
40	22	K	542 spiked	0.188	0.174	1
41	918	QC1	CRM [0.15/0.03]	0.149	0.033	1
42	0	B	Baseline	0.001	0.001	1

43	23	S	451	0.015	0.000	1
44	24	S	334	0.020	0.001	1
45	25	S	360	0.013	0.000	1
46	26	S	282	0.022	0.002	1
47	27	S	399	0.017	0.002	1
48	28	K	399 spiked	0.160	0.158	1
49	918	QC1	CRM [0.15/0.03]	0.149	0.035	1
50	0	B	Baseline	0.001	0.001	1
51	916	S	Blank	0.004	0.005	1
52	916	S	Blank	0.004	0.005	1
53	916	S	Blank	0.004	0.004	1
54	916	S	Blank	0.004	0.004	1
55	916	S	Blank	0.004	0.004	1
56	916	S	Blank	0.004	0.005	1
57	916	S	Blank	0.004	0.005	1
58	916	S	Blank	0.004	0.004	1
59	916	S	Blank	0.004	0.004	1
60	916	S	Blank	0.004	0.004	1
61	918	QC1	CRM [0.15/0.03]	0.149	0.034	1
62	0	B	Baseline	0.001	0.001	1
63	922	D	Drift	0.203	0.193	1
64	0	B	Final Base	0.001	0.001	1

Corrections

Channel	1	3
Baseline	done	done
Drift	none	none
Carryover	done	done
%	0.44	0.45

Spike Recovery


Blank spiked Added	0.150	0.149
Blank spiked Found	0.147	0.122
Recovery %	98.390	81.540
425 spiked Added	0.150	0.150
425 spiked Found	0.147	0.155
Recovery %	98.012	103.613
542 spiked Added	0.148	0.149
542 spiked Found	0.134	0.154
Recovery %	90.497	103.249
399 spiked Added	0.149	0.150
399 spiked Found	0.143	0.156
Recovery %	95.845	103.855

Calibration

Channel 1	NOx
-----------	-----


Type	Observed	Calculated	Target	% Error
1C	55361	0.201	0.200	0.68
2C	28536	0.098	0.100	-1.70
3C	15849	0.050	0.050	-0.87
4C	8559	0.022	0.025	-13.77
5C	5856	0.011	0.010	11.73
6C	4006	0.004	0.005	-18.68
7C	3985	0.004	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9994
 Carryover (%): 0.4356
 Calibration Equation: $y = bx + a$
 $y =$ Peak Height Digit. Units
 $x =$ Concentration Unit mg/L
 $a =$ 2.947547E+003
 $b =$ 2.603056E+005

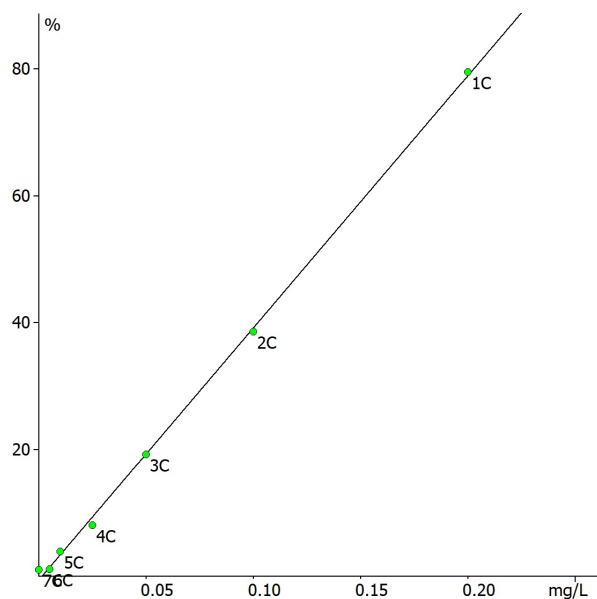

Type	Observed	Calculated	Target	% Error
1C	58108	0.199	0.200	-0.36
2C	31189	0.102	0.100	1.53
3C	17079	0.050	0.050	0.58
4C	9892	0.024	0.025	-3.24
5C	5934	0.010	0.010	-1.83
6C	4580	0.005	0.005	-1.99

Calibration Fit: Linear
 Correlation Coefficient: 0.9999
 Carryover (%): 0.4469
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = $3.230498E+003$
 b = $2.753821E+005$

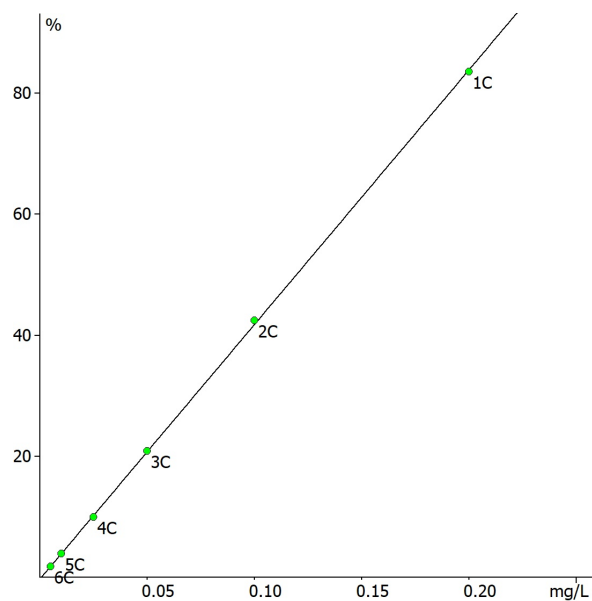

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero