

General

Name of Run	240716_Q5_NP1_DriftCorrected.run	Name of Analysis	NOx_P04.ANL
Date of Report	7/18/2024	System	AA500
Run Start	7/17/2024 1:15:10 PM	Run Stop	2:59:05 PM
Operator Run/Rec.	CLMR / CLMR	Software Version	AACE 8.02 SP2
Comment	Recalculate from Run240716_Q5_NP1.run		

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9991	1.0000	
AUFS/Gain	0.19/54	0.10/97	Data
Offset	246	4931	Field
Sensitivity	0.1308	0.0858	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.012	-0.011	0
1	921	P	Primer	0.243	0.200	1
2	921	D	Drift	0.208	0.200	1
3	921	C	0.200	0.199	0.200	1
4	901	C	0.100	0.100	0.100	1
5	902	C	0.050	0.053	0.050	1
6	903	C	0.025	0.029	0.025	1
7	904	C	0.010	0.009	0.010	1
8	0	X	Null-Cal.	0.006 N	0.010 N	1
9	906	C	Blank	-0.005	-0.001	1
10	0	B	Baseline	-0.012	-0.011	1
11	921	H1	High	0.231	0.200	1
12	905	L1	Low	0.003	0.010	1
13	905	L1	Low	0.005	0.010	1
14	906	S	Blank	0.000	-0.001	1
15	907	K	Blank spiked	0.166	0.146	1
16	908	QC1	CRM [0.15/0.03]	0.159	0.027	1
17	909	QC1	CRM 4/30	0.171	0.027	1
18	1	S	152	0.021	-0.002	1
19	2	S	347	0.030	0.020	1
20	3	S	100	0.021	0.000	1
21	4	S	204	0.021	0.000	1
22	5	S	87	0.021	-0.001	1
23	6	S	61	0.012	0.009	1
24	7	S	269	0.018	0.001	1
25	8	S	9	0.015	-0.003	1
26	9	S	243	0.013	0.004	1
27	10	S	165	0.022	-0.004	1
28	11	K	165 spiked	0.186	0.142	1
29	908	QC1	CRM [0.15/0.03]	0.155	0.027	1
30	909	QC1	CRM 4/30	0.166	0.027	1
31	0	B	Baseline	-0.012	-0.011	1
32	12	S	139	0.025	-0.003	1
33	13	S	35	0.014	0.003	1
34	14	S	22	0.006	-0.004	1
35	15	S	230	0.011	-0.004	1
36	16	S	48	0.025	-0.004	1
37	17	S	308	0.016	-0.004	1
38	18	S	191	0.011	-0.004	1
39	19	S	295	0.011	-0.003	1
40	20	S	282	0.012	-0.003	1
41	21	S	178	0.010	-0.004	1
42	22	K	178 spiked	0.152	0.141	1

43	908	QC1	CRM [0.15/0.03]	0.153	0.027	1
44	909	QC1	CRM 4/30	0.141	0.027	1
45	0	B	Baseline	-0.012	-0.011	1
46	23	S	74	0.009	0.002	1
47	24	S	113	0.013	-0.004	1
48	25	S	126	0.015	-0.005	1
49	26	S	321	0.016	-0.004	1
50	27	S	217	0.026	-0.005	1
51	28	S	256	0.020	-0.005	1
52	29	K	256 spiked	0.148	0.138	1
53	908	QC1	CRM [0.15/0.03]	0.167	0.026	1
54	909	QC1	CRM 4/30	0.163	0.026	1
55	0	B	Baseline	-0.012	-0.011	1
56	906	S	Blank	-0.006	-0.001	1
57	906	S	Blank	-0.005	0.000	1
58	906	S	Blank	-0.005	-0.002	1
59	906	S	Blank	-0.005	-0.001	1
60	906	S	Blank	-0.005	-0.002	1
61	906	S	Blank	-0.005	-0.002	1
62	906	S	Blank	-0.005	-0.002	1
63	906	S	Blank	-0.004	0.000	1
64	906	S	Blank	-0.004	-0.001	1
65	906	S	Blank	-0.003	0.000	1
66	908	QC1	CRM [0.15/0.03]	0.167	0.027	1
67	909	QC1	CRM 4/30	0.190	0.027	1
68	0	B	Baseline	-0.012	-0.011	1
69	921	D	Drift	0.208	0.200	1
70	0	B	Final Base	-0.012	-0.011	1

Corrections

Channel	1	3
Baseline	done	done
Drift	done	done
Carryover	done	done
%	9.54	0.38

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.166	0.147
Recovery %	110.416	97.863
165 spiked Added	0.149	0.150
165 spiked Found	0.164	0.146
Recovery %	109.732	97.125
178 spiked Added	0.150	0.150
178 spiked Found	0.142	0.145
Recovery %	95.159	96.849
256 spiked Added	0.149	0.150
256 spiked Found	0.128	0.144
Recovery %	85.776	95.636

Calibration

Channel 1 NOx

Type	Observed	Calculated	Target	% Error
1C	41269	0.199	0.200	-0.54
2C	23225	0.100	0.100	-0.38
3C	14796	0.053	0.050	6.45
4C	10407	0.029	0.025	16.29
5C	6777	0.009	0.010	-9.06
7C	4229	-0.005	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 0.9991
 Carryover (%): 9.5355
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 5.124674E+003
 b = 1.817012E+005

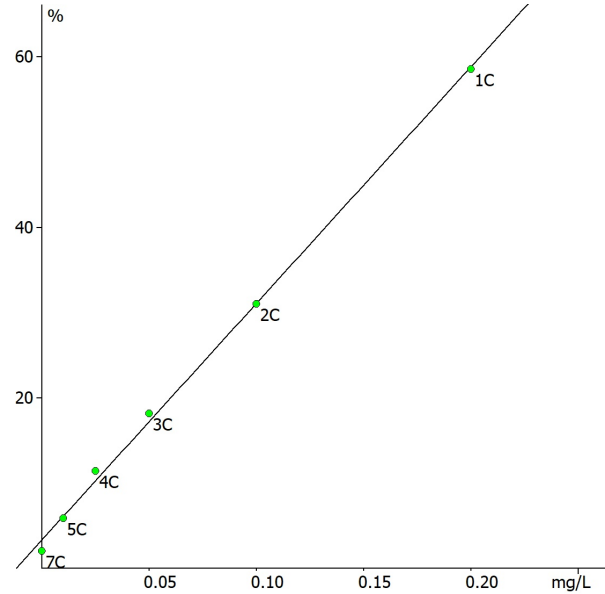

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	57837	0.200	0.200	-0.17
2C	32161	0.100	0.100	0.37
3C	19243	0.050	0.050	0.84
4C	12775	0.025	0.025	1.63
5C	8792	0.010	0.010	0.06
7C	5981	-0.001	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 1.0000
 Carryover (%): 0.3828
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 6.204299E+003
 b = 2.586057E+005

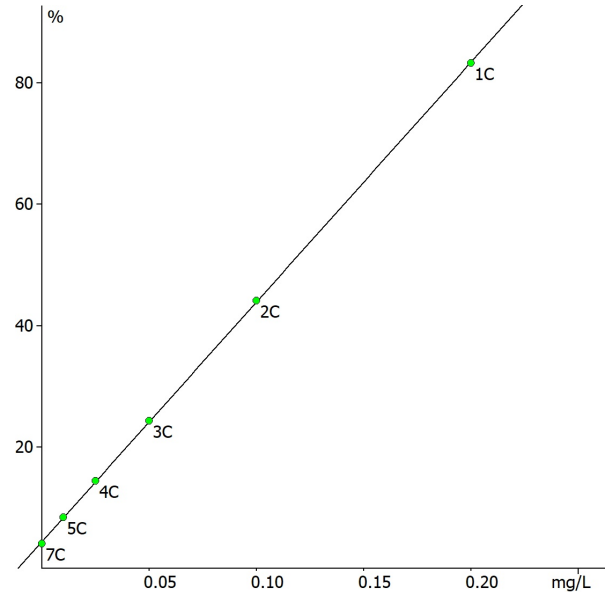

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero