

General

Name of Run	240801_Q5_TKN4_2C5CNull.run	Name of Analysis	TKN.ANL
Date of Report	8/1/2024	System	AA500
Run Start	8/1/2024 1:29:01 PM	Run Stop	4:12:51 PM
Operator Run/Rec.	CLMR / CLMR	Software Version	AACE 8.02 SP2
Comment	Recalculate from Run240801_Q5_TKN4R1.run te from Run240801_Q5_TKN4.run		

Parameters

Channel	1	
Method	TKN	Dilution
Unit	mg/L	
Calib. Fit	Linear	
Corr.Coeff.(r)	0.9990	
AUFS/Gain	0.16/63	Data
Offset	-800	Field
Sensitivity	0.1288	'#2
Sample Limit 1	-	
Sample Limit 2	-	

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.432	0	
1	921	P	Primer	2.224	1	
2	921	D	Drift	2.225	1	
3	921	C	2.22	2.227	1	
4	0	X		1.277 N	1	
5	902	C	1.11	1.068	1	
6	903	C	0.55	0.606	1	
7	0	X		0.635 N	1	
8	905	C	Blank	-0.022	1	
9	0	B	Baseline	-0.432	1	
10	921	H1	High	2.227	1	
11	904	L1	Low	0.635	1	
12	904	L1	Low	0.635	1	
13	905	S	Blank	-0.025	1	
14	906	K	Blank Spike	0.948	1	
15	907	QC1	CRM [1.5 mg/L]	1.483	1	
16	1	S	10	0.086	1	
17	2	S	75	0.399	1	
18	3	S	114	0.212	1	
19	4	S	62	0.115	1	
20	5	S	88	0.140	1	
21	6	S	205	0.086	1	
22	7	S	270	0.962	1	
23	8	S	101	0.158	1	
24	9	S	127	0.080	1	
25	10	S	439	0.034	1	
26	11	K	439 spiked	1.217	1	
27	907	QC1	CRM [1.5 mg/L]	1.532	1	
28	0	B	Baseline	-0.432	1	
29	12	S	153	0.125	1	
30	13	S	283	0.071	1	
31	14	S	322	0.554	1	
32	15	S	348	0.154	1	
33	16	S	166	0.295	1	
34	17	S	140	0.409	1	
35	18	S	23	0.169	1	
36	19	S	36	0.374	1	
37	20	S	179	0.341	1	
38	21	S	569	0.116	1	
39	22	K	569 spiked	1.142	1	
40	907	QC1	CRM [1.5 mg/L]	1.481	1	
41	0	B	Baseline	-0.432	1	
42	23	S	231	0.404	1	

43	24	S	218	0.117	1
44	25	S	257	-0.054	1
45	26	S	309	0.652	1
46	27	S	192	0.383	1
47	28	S	296	0.131	1
48	29	S	244	0.306	1
49	30	S	49	0.187	1
50	31	S	582	-0.012	1
51	32	K	582 spiked	1.135	1
52	907	QC1	CRM [1.5 mg/L]	1.483	1
53	0	B	Baseline	-0.432	1
54	33	S	517	0.142	1
55	34	S	452	-0.012	1
56	35	S	400	0.017	1
57	36	S	413	0.088	1
58	37	S	374	-0.018	1
59	38	S	543	0.076	1
60	39	S	504	-0.014	1
61	40	S	608	-0.017	1
62	41	S	426	0.006	1
63	42	S	361	0.032	1
64	43	K	361 spiked	1.160	1
65	907	QC1	CRM [1.5 mg/L]	1.477	1
66	0	B	Baseline	-0.432	1
67	44	S	621	0.020	1
68	45	S	478	0.000	1
69	46	S	335	0.188	1
70	47	S	530	0.305	1
71	48	S	556	0.035	1
72	49	S	387	-0.008	1
73	50	S	491	0.001	1
74	51	S	335	1.036	1
75	52	S	465	0.120	1
76	53	K	465 spiked	1.214	1
77	907	QC1	CRM [1.5 mg/L]	1.466	1
78	0	B	Baseline	-0.432	1
79	54	S	348 spiked IW	0.977	1
80	55	S	23 spiked IW	0.998	1
81	56	S	270 spiked IW	-0.015	1
82	57	S	296 spiked CW	0.870	1
83	58	S	244 spiked CW	0.902	1
84	59	S	49 spiked CW	0.832	1
85	907	QC1	CRM [1.5 mg/L]	1.474	1
86	0	B	Baseline	-0.432	1
87	905	S	Blank	-0.028	1
88	905	S	Blank	-0.032	1
89	905	S	Blank	-0.030	1
90	905	S	Blank	-0.029	1
91	905	S	Blank	-0.022	1
92	905	S	Blank	-0.035	1
93	905	S	Blank	-0.026	1
94	905	S	Blank	-0.031	1
95	905	S	Blank	-0.034	1
96	905	S	Blank	-0.034	1
97	906	K	Blank spiked	0.934	1
98	907	QC1	CRM [1.5 mg/L]	1.466	1
99	0	B	Baseline	-0.432	1
100	921	D	Drift	2.225	1
101	0	B	Final Base	-0.432	1

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	done
%	0.09

Spike Recovery

Blank Spike Added	1.001
Blank Spike Found	0.974
Recovery %	97.286
439 spiked Added	1.001
439 spiked Found	1.183


Recovery %	118.216
569 spiked Added	1.000
569 spiked Found	1.026
Recovery %	102.642
582 spiked Added	1.001
582 spiked Found	1.147
Recovery %	114.552
361 spiked Added	1.001
361 spiked Found	1.128
Recovery %	112.737
465 spiked Added	1.000
465 spiked Found	1.094
Recovery %	109.404
Blank spiked Added	1.001
Blank spiked Found	0.967
Recovery %	96.618

Calibration

Channel 1 TKN

Type	Observed	Calculated	Target	% Error
1C	57109	2.227	2.220	0.31
3C	33904	1.068	1.110	-3.74
4C	24644	0.606	0.550	10.23
6C	12066	-0.022	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9990
 Carryover (%): 0.0911
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 1.249830E+004
 b = 2.003346E+004

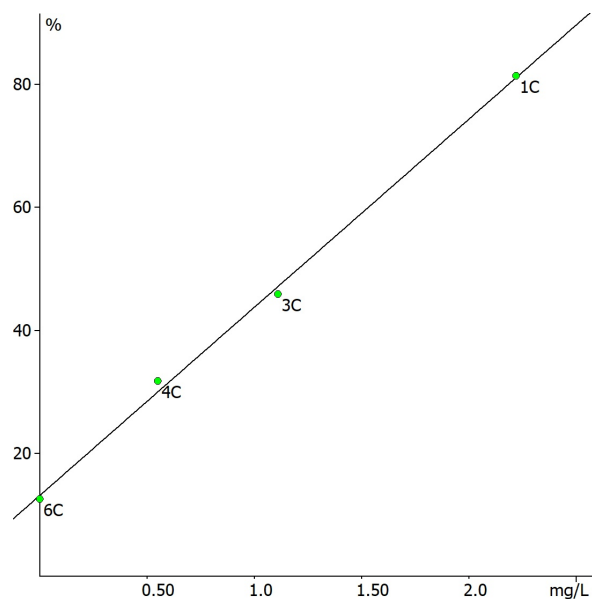

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero