

General

Name of Run	240805_Q5_NP2R1.RUN	Name of Analysis	NOx_P04.ANL
Date of Report	8/5/2024	System	AA500
Run Start	8/5/2024 1:55:33 PM	Run Stop	3:35:13 PM
Operator Run/Rec. CLMR	/	Software Version	AACE 8.02 SP2
Comment			

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9990	0.9998	
AUFS/Gain	0.22/46	0.02/603	Data
Offset	239	5051	Field
Sensitivity	0.1508	0.0140	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.005	-0.008	0
1	921	P	Primer	0.200	0.198	1
2	921	D	Drift	0.199	0.199	1
3	921	C	0.200	0.198 M	0.200	1
4	901	C	0.100	0.105	0.100	1
5	902	C	0.050	0.046	0.053	1
6	903	C	0.025	0.024	0.024	1
7	904	C	0.010	0.012	0.010	1
8	905	C	0.005	0.006	0.004	1
9	906	C	Blank	-0.002	0.000	1
10	0	B	Baseline	-0.005	-0.008	1
11	921	H1	High	0.195	0.198	1
12	905	L1	Low	0.004	0.002	1
13	905	L1	Low	0.006	0.002	1
14	906	S	Blank	-0.001	0.001	1
15	907	K	Blank spiked	0.123	0.142	1
16	908	QC1	CRM 7/15 [0.15/0.03]	0.143	0.027	1
17	909	QC1	CRM 4/30	0.152	0.027	1
18	1	S	412	0.010	0.002	1
19	2	S	451	0.011	0.003	1
20	3	S	477	0.011	0.003	1
21	4	S	529	0.020	0.003	1
22	5	S	360	0.004	0.002	1
23	6	S	490	0.011	0.002	1
24	7	S	464	0.011	0.004	1
25	8	S	542	-0.003	0.001	1
26	9	S	399	0.011	0.003	1
27	10	S	425	0.007	0.001	1
28	11	K	425 spiked	0.151	0.150	1
29	908	QC1	CRM 7/15 [0.15/0.03]	0.152	0.029	1
30	909	QC1	CRM 4/30	0.146	0.029	1
31	0	B	Baseline	-0.005	-0.008	1
32	12	S	386	0.011	0.003	1
33	13	S	581	0.011	0.006	1
34	14	S	334	0.010	0.003	1
35	15	S	373	0.008	0.004	1
36	16	S	516	0.005	0.002	1
37	17	S	568	0.017	0.003	1
38	18	S	607	0.009	0.005	1
39	19	S	620	0.010	0.002	1
40	20	S	555	0.008	0.003	1
41	21	S	503	0.008	0.004	1
42	22	K	503 spiked	0.129	0.149 M	1

43	23	S	438	0.012	0.005	1
44	24	K	438 spiked	0.221	0.150	1
45	908	QC1	CRM 7/15 [0.15/0.03]	0.173	0.030	1
46	909	QC1	CRM 4/30	0.161	0.029	1
47	0	B	Baseline	-0.005	-0.008	1
48	906	S	Blank	-0.002	-0.002	1
49	906	S	Blank	-0.003	0.000	1
50	906	S	Blank	-0.003	-0.001	1
51	906	S	Blank	-0.003	0.000	1
52	906	S	Blank	-0.004	-0.002	1
53	906	S	Blank	-0.004	-0.002	1
54	906	S	Blank	-0.004	-0.002	1
55	906	S	Blank	-0.004	0.000	1
56	906	S	Blank	-0.005	-0.002	1
57	906	S	Blank	-0.002	-0.002	1
58	907	K	Blank spiked	0.156	0.143	1
59	908	QC1	CRM 7/15 [0.15/0.03]	0.174	0.029	1
60	909	QC1	CRM 4/30	0.176	0.029	1
61	0	B	Baseline	-0.005	-0.008	1
62	25	S	607 spiked	0.182	0.148	1
63	26	S	620 spiked	0.174	0.152	1
64	27	S	542	0.029	0.040	1
65	0	B	Baseline	-0.005	-0.008	1
66	921	D	Drift	0.199	0.199	1
67	0	B	Final Base	-0.005	-0.008	1

Corrections

Channel	1	3
Baseline	done	done
Drift	done	done
Carryover	done	none
%	7.61	-0.02

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.124	0.141
Recovery %	82.948	93.770
425 spiked Added	0.150	0.150
425 spiked Found	0.143	0.149
Recovery %	95.727	99.571
503 spiked Added	0.150	0.150
503 spiked Found	0.121	0.145
Recovery %	81.036	96.767
438 spiked Added	0.150	0.150
438 spiked Found	0.209	0.144
Recovery %	139.908	96.334
Blank spiked Added	0.150	0.150
Blank spiked Found	0.158	0.144
Recovery %	105.407	96.195

Calibration

Channel 1 NOx

Type	Observed	Calculated	Target	% Error
1C	48144	0.198	0.200	-0.81
2C	27590	0.105	0.100	5.47
3C	14351	0.046	0.050	-8.76
4C	9498	0.024	0.025	-5.29
5C	6936	0.012	0.010	20.97
6C	5692	0.006	0.005	29.46
7C	3880	-0.002	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 0.9990
 Carryover (%): 7.6085
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 4.260158E+003
 b = 2.212048E+005

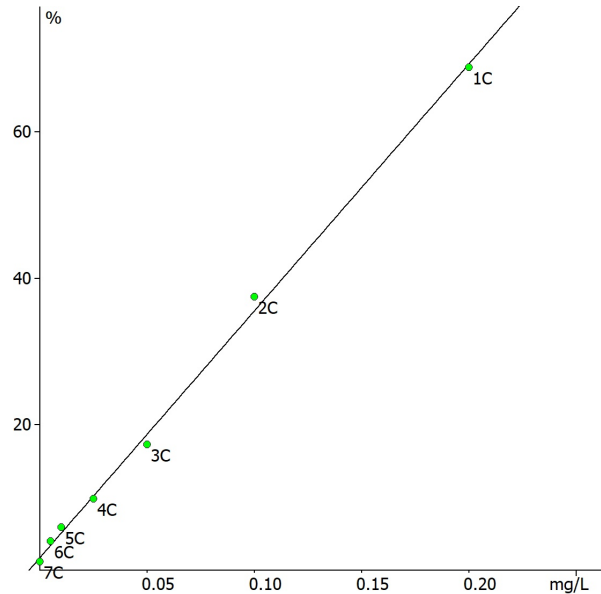

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	59140	0.200	0.200	-0.24
2C	32454	0.100	0.100	-0.14
3C	19829	0.053	0.050	5.43
4C	12162	0.024	0.025	-3.65
5C	8488	0.010	0.010	3.67
6C	6732	0.004	0.005	-23.81
7C	5616	0.000	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9998
 Carryover (%): no
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 5.711836E+003
 b = 2.677901E+005

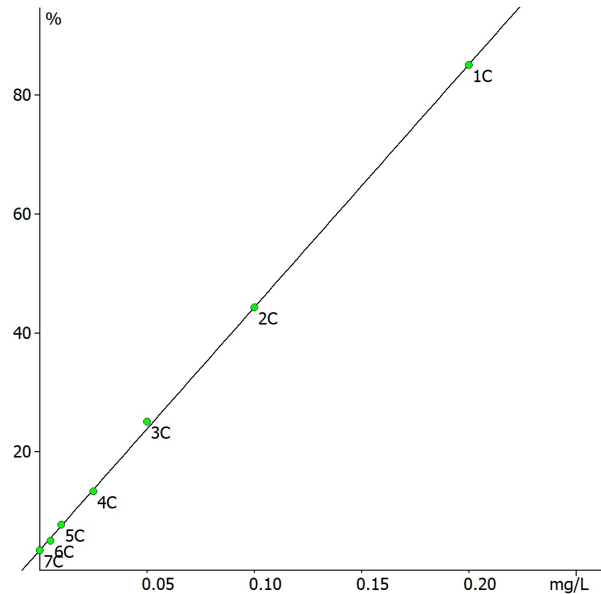

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero