

General

Name of Run	240918_Q6_TKN1.run	Name of Analysis	TKN.ANL
Date of Report	9/19/2024	System	AA500
Run Start	9/18/2024 3:46:12 PM	Run Stop	5:29:54 PM
Operator Run/Rec.	CLMR_CVMD	Software Version	AACE 8.02 SP2
Comment	1st Batch		

Parameters

Channel	1		
Method	TKN	Dilution	
Unit	mg/L		
Calib. Fit	Linear		
Corr.Coeff.(r)	0.9985		
AUFS/Gain	0.16/62	Data	
Offset	-1069	Field	
Sensitivity	0.1372	'#2	
Sample Limit 1	-		
Sample Limit 2	-		

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.325	0	
1	921	P	Primer	2.278	1	
2	921	D	Drift	2.280	1	
3	921	C	2.22	2.279	1	
4	901	C	1.66	1.582	1	
5	902	C	1.11	1.119	1	
6	903	C	0.55	0.518	1	
7	904	C	0.22	0.241	1	
8	905	C	Blank	0.020	1	
9	0	B	Baseline	-0.325	1	
10	921	H1	High	2.281	1	
11	904	L1	Low	0.237	1	
12	904	L1	Low	0.237	1	
13	905	S	Blank	0.019	1	
14	906	K	Blank Spike	0.969	1	
15	907	QC1	CRM [1.5 mg/L]	1.816	1	
16	1	S	49	0.080	1	
17	2	S	88	0.104	1	
18	3	S	140	0.128	1	
19	4	S	127	0.092	1	
20	5	S	23	0.114	1	
21	6	S	62	0.078	1	
22	7	S	244	0.072	1	
23	8	S	75	0.071	1	
24	9	S	101	0.179	1	
25	10	S	179	0.051	1	
26	11	K	179 spiked	1.039	1	
27	907	QC1	CRM [1.5 mg/L]	1.815	1	
28	0	B	Baseline	-0.325	1	
29	12	S	218	0.047	1	
30	13	S	205	0.071	1	
31	14	S	166	0.102	1	
32	15	S	296	0.054	1	
33	16	S	270	0.119	1	
34	17	S	322	0.087	1	
35	18	S	114	0.114	1	
36	19	S	153	0.088	1	
37	20	S	283	0.155	1	
38	21	S	36	0.075	1	
39	22	K	36 spiked	1.070	1	
40	907	QC1	CRM [1.5 mg/L]	1.818	1	
41	0	B	Baseline	-0.325	1	
42	23	S	335	0.052	1	

43	24	S	192	0.074	1
44	25	S	10	0.080	1
45	26	K	10 spiked	1.042	1
46	907	QC1	CRM [1.5 mg/L]	1.817	1
47	0	B	Baseline	-0.325	1
48	905	S	Blank	0.023	1
49	905	S	Blank	0.021	1
50	905	S	Blank	0.018	1
51	905	S	Blank	0.021	1
52	905	S	Blank	0.018	1
53	905	S	Blank	0.017	1
54	905	S	Blank	0.016	1
55	905	S	Blank	0.018	1
56	905	S	Blank	0.018	1
57	905	S	Blank	0.016	1
58	906	K	Blank spiked	0.963	1
59	907	QC1	CRM [1.5 mg/L]	1.813	1
60	0	B	Baseline	-0.325	1
61	921	D	Drift	2.280	1
62	0	B	Final Base	-0.325	1

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	done
%	0.17

Spike Recovery


Blank Spike Added	1.001
Blank Spike Found	0.950
Recovery %	94.936
179 spiked Added	1.000
179 spiked Found	0.988
Recovery %	98.715
36 spiked Added	1.000
36 spiked Found	0.994
Recovery %	99.400
10 spiked Added	1.000
10 spiked Found	0.962
Recovery %	96.168
Blank spiked Added	1.001
Blank spiked Found	0.947
Recovery %	94.668

Calibration

Channel 1	TKN
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	58921	2.279	2.220	2.67
2C	43993	1.582	1.660	-4.68
3C	34081	1.119	1.110	0.85
4C	21194	0.518	0.550	-5.87
5C	15266	0.241	0.220	9.50
6C	10543	0.020	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9985
 Carryover (%): 0.1718
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 1.010700E+004
 b = 2.141597E+004

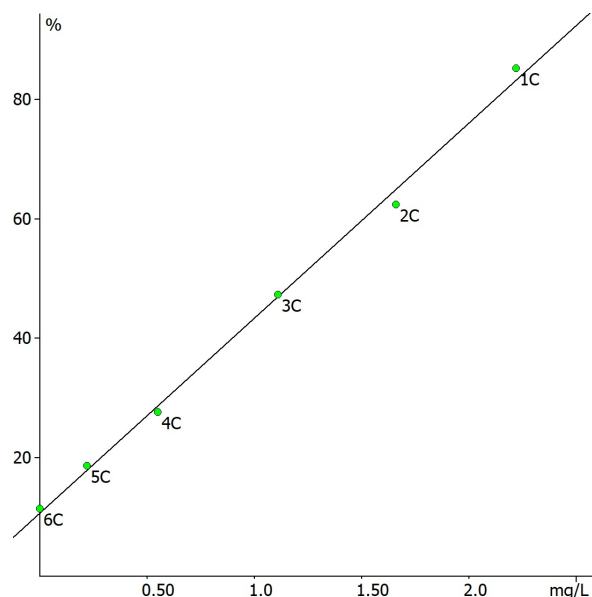

Legend

* ... Sample offscale
+ ... Result higher than sample limit
- ... Result lower than sample limit
P ... Standard passed
F ... Standard failed
N ... Value not calculated or not used
R ... Resample after offscale
M ... Peak marker moved manually
d ... Pre-Diluted sample
D ... Post-Diluted sample
H ... Dual calibration used
Z ... Report negative values as zero