

General

Name of Run	241002_Q6_TKN6.run	Name of Analysis	TKN.ANL
Date of Report	10/2/2024	System	AA500
Run Start	10/2/2024 3:10:11 PM	Run Stop	4:40:17 PM
Operator Run/Rec.	CLMR	Software Version	AACE 8.02 SP2
Comment	2nd Batch		
	CalCurve_TKN1		
	240918_Q6_TKN1.run		

Parameters

Channel	1		
Method	TKN	Dilution	
Unit	mg/L		
Calib. Fit	Linear		
Corr.Coeff.(r)	0.9985		
AUFS/Gain	0.14/73	Data	
Offset	-1395	Field	
Sensitivity	0.1156	'#2	
Sample Limit 1	-		
Sample Limit 2	-		

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.325	0	
1	921	P	Primer	2.257	1	
2	921	D	Drift	2.257	1	
3	905	C	Blank	-0.107	1	
4	0	B	Baseline	-0.325	1	
5	921	H1	High	2.253	1	
6	904	L1	Low	0.580	1	
7	904	L1	Low	0.580	1	
8	905	S	Blank	-0.101	1	
9	906	K	Blank Spike	0.962	1	
10	908	QC1	CRM 9/24 [1.5 mg/L]	1.635	1	
11	0	B	Baseline	-0.325	1	
12	1	S	556	-0.153	1	
13	2	S	452	0.009	1	
14	3	S	426	-0.173	1	
15	4	S	504	-0.144	1	
16	5	S	517	-0.154	1	
17	6	S	530	-0.197	1	
18	7	S	595	-0.061	1	
19	8	S	569	-0.141	1	
20	9	S	400	-0.142	1	
21	10	S	413	-0.145	1	
22	11	S	543	-0.145	1	
23	12	S	387	-0.164	1	
24	13	K	387 spiked	1.110	1	
25	908	QC1	CRM 9/24 [1.5 mg/L]	1.614	1	
26	0	B	Baseline	-0.325	1	
27	14	S	309	-0.209	1	
28	15	S	374	-0.163	1	
29	16	S	465	-0.206	1	
30	17	S	348	-0.225	1	
31	18	S	361	-0.167	1	
32	19	S	491	-0.136	1	
33	20	S	231	-0.172	1	
34	21	S	439	-0.235	1	
35	22	S	478	-0.221	1	
36	23	S	608	0.104	1	
37	24	S	582	-0.246	1	
38	25	K	582 spiked	0.680	1	
39	908	QC1	CRM 9/24 [1.5 mg/L]	1.606	1	
40	0	B	Baseline	-0.325	1	
41	921	D	Drift	2.257	1	

42	905	S	Blank	-0.106	1
43	905	S	Blank	-0.116	1
44	905	S	Blank	-0.127	1
45	905	S	Blank	-0.133	1
46	905	S	Blank	-0.141	1
47	906	K	Blank spiked	0.967	1
48	908	QC1	CRM 9/24 [1.5 mg/L]	1.594	1
49	0	B	Baseline	-0.325	1
50	921	D	Drift	2.257	1
51	0	B	Final Base	-0.325	1

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	done
%	0.13

Spike Recovery


Blank Spike Added	1.002
Blank Spike Found	1.063
Recovery %	106.100
387 spiked Added	1.002
387 spiked Found	1.274
Recovery %	127.134
582 spiked Added	1.003
582 spiked Found	0.926
Recovery %	92.362
Blank spiked Added	1.002
Blank spiked Found	1.108
Recovery %	110.581

Calibration

Channel 1	TKN
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	58921	2.279	2.220	2.67
2C	43993	1.582	1.660	-4.68
3C	34081	1.119	1.110	0.85
4C	21194	0.518	0.550	-5.87
5C	15266	0.241	0.220	9.50
6C	10543	0.020	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9985
 Carryover (%): 0.1315
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 1.010700E+004
 b = 2.141597E+004

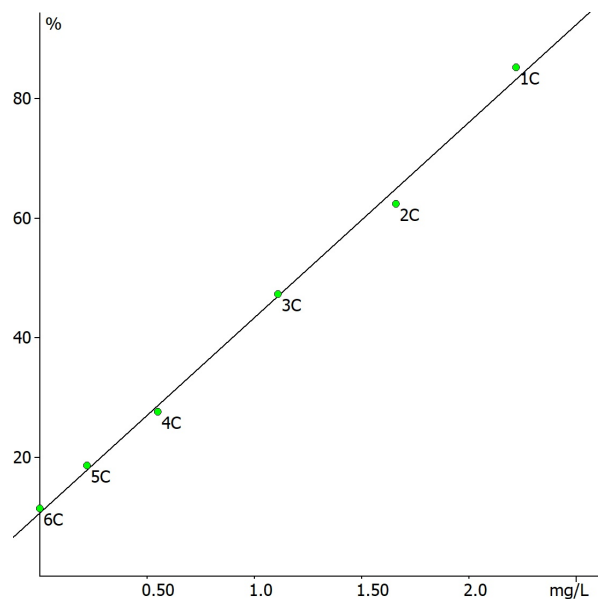

Legend

* ... Sample offscale
+ ... Result higher than sample limit
- ... Result lower than sample limit
P ... Standard passed
F ... Standard failed
N ... Value not calculated or not used
R ... Resample after offscale
M ... Peak marker moved manually
d ... Pre-Diluted sample
D ... Post-Diluted sample
H ... Dual calibration used
Z ... Report negative values as zero