

General

Name of Run	241212_Q7_NP1.run	Name of Analysis	NOx_P04.ANL
Date of Report	12/12/2024	System	AA500
Run Start	12/12/2024 2:22:45 PM	Run Stop	4:07:37 PM
Operator Run/Rec.	CLMR&CVMD	Software Version	AACE 8.02 SP2
Comment			

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9999	0.9999	
AUFS/Gain	0.45/22	0.10/101	Data
Offset	158	4879	Field
Sensitivity	0.3718	0.0832	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.003	-0.006	0
1	921	P	Primer	0.199	0.198	1
2	921	D	Drift	0.198	0.199	1
3	921	C	0.200	0.200	0.199	1
4	901	C	0.100	0.099	0.101	1
5	902	C	0.050	0.050	0.051	1
6	903	C	0.025	0.026	0.024	1
7	904	C	0.010	0.010	0.010	1
8	905	C	0.005	0.006	0.005	1
9	906	C	Blank	-0.001	-0.001	1
10	0	B	Baseline	-0.003	-0.006	1
11	921	H1	High	0.200	0.198	1
12	905	L1	Low	0.006	0.005	1
13	905	L1	Low	0.006	0.005	1
14	906	S	Blank	-0.001	-0.001	1
15	907	K	Blank spiked	0.144	0.142	1
16	908	QC1	CRM 7/15 [0.15/0.03]	0.153	0.030	1
17	909	QC1	CRM 4/30 [0.15/0.03]	0.153	0.030	1
18	1	S	74	0.008	0.001	1
19	2	S	100	0.013	0.002	1
20	3	S	48	0.005	0.001	1
21	4	S	191	0.012	0.001	1
22	5	S	35	0.009	0.001	1
23	6	S	152	0.009	0.000	1
24	7	S	139	0.008	0.000	1
25	8	S	61	0.004	0.000	1
26	9	S	217	0.013	0.000	1
27	10	S	9	0.017	0.001	1
28	11	K	9 spiked	0.165	0.148	1
29	908	QC1	CRM 7/15 [0.15/0.03]	0.154	0.030	1
30	909	QC1	CRM 4/30 [0.15/0.03]	0.154	0.030	1
31	0	B	Baseline	-0.003	-0.006	1
32	12	S	126	0.006	0.001	1
33	13	S	22	0.007	0.001	1
34	14	S	87	0.005	0.000	1
35	15	S	256	0.008	0.000	1
36	16	S	178	0.010	0.000	1
37	17	S	165	0.004	0.001	1
38	18	S	204	0.009	0.000	1
39	19	S	269	0.007	0.001	1
40	20	S	113	0.003	0.001	1
41	21	S	386	0.004	0.001	1
42	22	K	386 spiked	0.148	0.144	1

43	908	QC1	CRM 7/15 [0.15/0.03]	0.151	0.031	1
44	909	QC1	CRM 4/30 [0.15/0.03]	0.152	0.030	1
45	0	B	Baseline	-0.003	-0.006	1
46	23	S	295	0.003	0.000	1
47	24	S	308	0.007	-0.001	1
48	25	S	464	0.007	0.002	1
49	26	S	334	0.007	0.001	1
50	27	S	282	0.006	0.001	1
51	28	S	477	0.008	0.001	1
52	29	S	399	0.005	0.000	1
53	30	K	399 spiked	0.151	0.145	1
54	908	QC1	CRM 7/15 [0.15/0.03]	0.154	0.034	1
55	909	QC1	CRM 4/30 [0.15/0.03]	0.151	0.031	1
56	0	B	Baseline	-0.003	-0.006	1
57	906	S	Blank	-0.001	-0.001	1
58	906	S	Blank	-0.001	-0.002	1
59	906	S	Blank	-0.002	-0.002	1
60	906	S	Blank	-0.002	-0.002	1
61	906	S	Blank	-0.002	-0.003	1
62	906	S	Blank	-0.002	-0.003	1
63	906	S	Blank	-0.002	-0.003	1
64	906	S	Blank	-0.001	-0.003	1
65	906	S	Blank	-0.001	-0.003	1
66	906	S	Blank	-0.001	-0.003	1
67	907	K	Blank spiked	0.141	0.143	1
68	908	QC1	CRM 7/15 [0.15/0.03]	0.150	0.031	1
69	909	QC1	CRM 4/30 [0.15/0.03]	0.150	0.031	1
70	23	S	321	0.003	0.001	1
71	0	B	Baseline	-0.003	-0.006	1
72	921	D	Drift	0.198	0.199	1
73	0	B	Final Base	-0.003	-0.006	1

Corrections

Channel	1	3
Baseline	done	done
Drift	done	done
Carryover	done	done
%	0.61	0.03

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.145	0.143
Recovery %	96.674	95.623
9 spiked Added	0.149	0.150
9 spiked Found	0.147	0.147
Recovery %	98.625	97.730
386 spiked Added	0.150	0.150
386 spiked Found	0.144	0.144
Recovery %	96.263	95.803
399 spiked Added	0.150	0.150
399 spiked Found	0.147	0.144
Recovery %	97.936	96.186
Blank spiked Added	0.150	0.150
Blank spiked Found	0.142	0.146
Recovery %	94.778	97.181

Calibration

Channel 1 NOx

Type	Observed	Calculated	Target	% Error
1C	57131	0.200	0.200	0.08
2C	30420	0.099	0.100	-0.62
3C	17382	0.050	0.050	0.37
4C	10928	0.026	0.025	3.33
5C	6604	0.010	0.010	-4.83
6C	5723	0.006	0.005	23.87
7C	3748	-0.001	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 0.9999
 Carryover (%): 0.6097
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 4.081494E+003
 b = 2.650398E+005

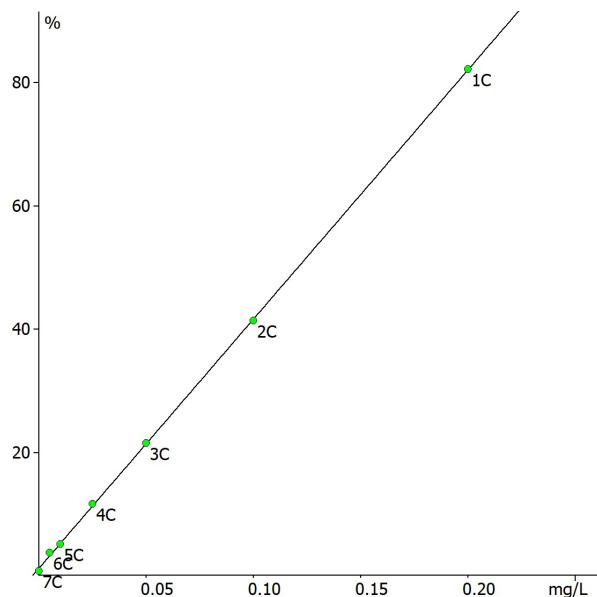

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58650	0.199	0.200	-0.40
2C	32109	0.101	0.100	1.29
3C	18441	0.051	0.050	1.72
4C	11261	0.024	0.025	-2.52
5C	7447	0.010	0.010	2.99
6C	5985	0.005	0.005	-1.89
7C	4403	-0.001	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9999
 Carryover (%): 0.0287
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 4.655349E+003
 b = 2.710476E+005

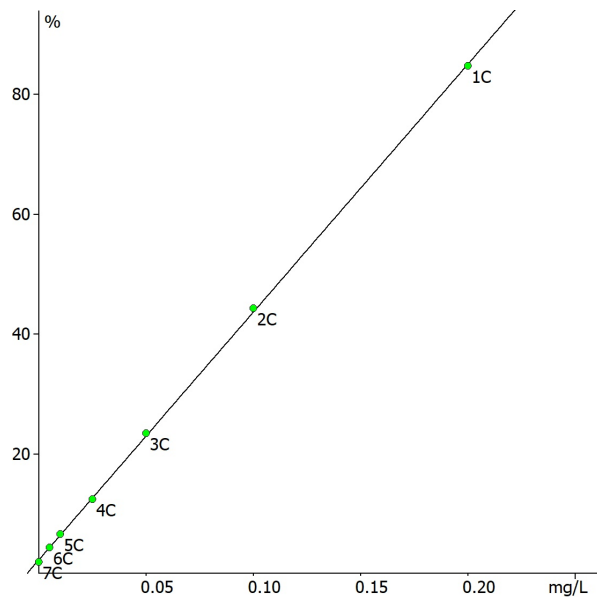

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero