

General

Name of Run	250203_Q7_NP4R1_PhosphateCurveTAC.run	Name of Analysis	NOx_P04.ANL
Date of Report	2/7/2025	System	AA500
Run Start	2/3/2025 3:58:40 PM	Run Stop	5:17:08 PM
Operator Run/Rec.	CLMR / tac	Software Version	AACE 8.02 SP2
Comment	Recalculate from Run250203_Q7_NP4.run		

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9997	1.0000	
AUFS/Gain	0.38/26	0.14/74	Data
Offset	149	4775	Field
Sensitivity	0.3648	0.1125	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.003	-0.003	0
1	921	P	Primer	0.189	0.199	1
2	921	D	Drift	0.192 *	0.200	1
3	921	C	0.200	0.198 *	0.200	1
4	0	XC	Cal. (1)	0.103	0.122 N	1
5	0	XC	Cal. (1)	0.050	0.081 N	1
6	0	XC	Cal. (1)	0.026	0.061 N	1
7	0	XC	Cal. (1)	0.009	0.048 N	1
8	0	XC	Cal. (1)	0.004	0.043 N	1
9	906	C	Blank	-0.001	0.000	1
10	0	B	Baseline	-0.003	-0.003	1
11	921	H1	High	0.197 *	0.199	1
12	905	L1	Low	0.005	0.043	1
13	905	L1	Low	0.005	0.043	1
14	906	S	Blank	-0.001	0.000	1
15	907	K	Blank spiked	0.145	0.117	1
16	910	QC1	CRM 2/3 [0.15/0.03]	0.148	0.024	1
17	1	S	633	0.005	0.011	1
18	2	S	360	0.008	0.004	1
19	3	S	542	0.007	-0.003	1
20	4	S	425	0.008	0.010	1
21	5	S	672	0.007	0.007	1
22	6	S	373	0.009	0.003	1
23	7	S	620	0.008	0.009	1
24	8	S	451	0.006	0.003	1
25	9	S	347	0.008	0.006	1
26	10	S	607	0.013	0.014	1
27	11	K	607 spiked	0.157	0.143	1
28	910	QC1	CRM 2/3 [0.15/0.03]	0.145	0.023	1
29	0	B	Baseline	-0.003	-0.003	1
30	12	S	659	0.009	0.008	1
31	13	S	490	0.007	0.003	1
32	14	S	529	0.010	0.003	1
33	15	S	685	0.006	0.006	1
34	16	S	646	0.016	0.020	1
35	17	K	646 spiked	0.164	0.158	1
36	910	QC1	CRM 2/3 [0.15/0.03]	0.142	0.023	1
37	0	B	Baseline	-0.003	-0.003	1
38	906	S	Blank	-0.002	-0.002	1
39	906	S	Blank	-0.002	-0.001	1
40	906	S	Blank	-0.002	-0.002	1
41	906	S	Blank	-0.002	-0.002	1
42	906	S	Blank	-0.002	-0.001	1

43	906	S	Blank	-0.002	-0.001	1
44	906	S	Blank	-0.002	-0.001	1
45	906	S	Blank	-0.002	-0.001	1
46	906	S	Blank	-0.002	-0.001	1
47	906	S	Blank	-0.001	0.000	1
48	907	K	Blank spiked	0.143	0.116	1
49	910	QC1	CRM 2/3 [0.15/0.03]	0.146	0.024	1
50	0	B	Baseline	-0.003	-0.003	1
51	921	D	Drift	0.194 *	0.200	1
52	0	B	Final Base	-0.003	-0.003	1

Corrections

Channel	1	3
Baseline	done	done
Drift	none	done
Carryover	none	done
%	0.00	0.38

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.147	0.117
Recovery %	97.782	78.115
607 spiked Added	0.150	0.150
607 spiked Found	0.144	0.129
Recovery %	96.491	86.250
646 spiked Added	0.150	0.149
646 spiked Found	0.148	0.138
Recovery %	98.752	92.080
Blank spiked Added	0.150	0.150
Blank spiked Found	0.145	0.116
Recovery %	96.384	77.584

Calibration

Channel 1	NOx
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	68479	0.198	0.200	-0.83
2C	37608	0.103	0.100	2.96
3C	20567	0.050	0.050	0.62
4C	12790	0.026	0.025	5.13
5C	7222	0.009	0.010	-9.20
6C	5713	0.004	0.005	-11.64
7C	3831	-0.001	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 0.9997
 Carryover (%): 0.0000
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 4.282998E+003
 b = 3.236642E+005

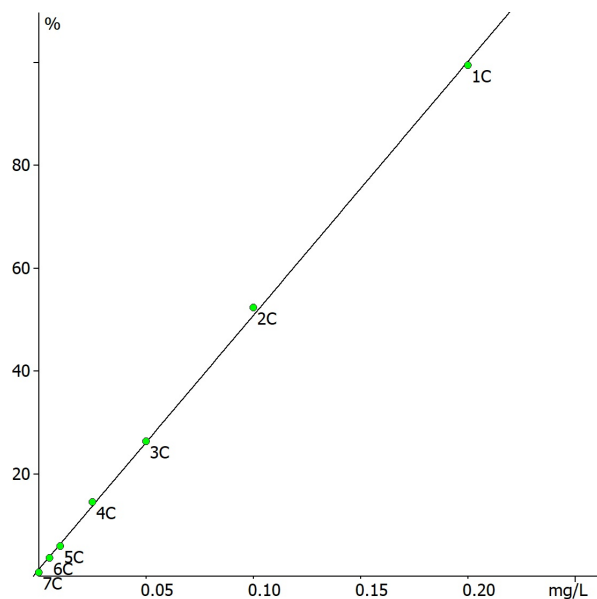

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58865	0.200	0.200	0.00
7C	4771	0.000	0.000	-

Calibration Fit: Linear
Correlation Coefficient: 1.0000
Carryover (%): 0.3825
Calibration Equation: $y = bx + a$
y = Peak Height Digit. Units
x = Concentration Unit mg/L
a = 4.771000E+003
b = 2.704700E+005

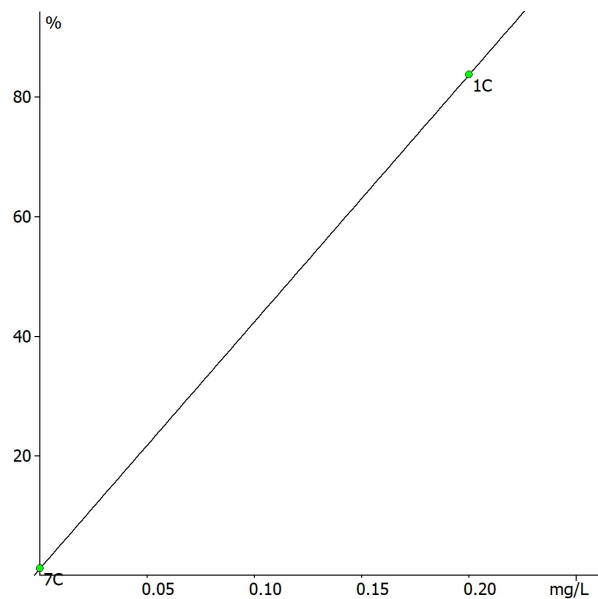

Legend

* ... Sample offscale
+ ... Result higher than sample limit
- ... Result lower than sample limit
P ... Standard passed
F ... Standard failed
N ... Value not calculated or not used
R ... Resample after offscale
M ... Peak marker moved manually
d ... Pre-Diluted sample
D ... Post-Diluted sample
H ... Dual calibration used
Z ... Report negative values as zero