

General

Name of Run	250207_Q7_TKN6.run	Name of Analysis	TKN.ANL
Date of Report	2/7/2025	System	AA500
Run Start	2/7/2025 1:28:50 PM	Run Stop	3:43:03 PM
Operator Run/Rec.	CLMR	Software Version	AACE 8.02 SP2
Comment	2nd Batch		

241206_Q7_TKN1.run

Parameters

Channel	1		
Method	TKN	Dilution	
Unit	mg/L		
Calib. Fit	Linear		
Corr.Coeff.(r)	0.9979		
AUFS/Gain	0.15/65	Data	
Offset	-1534	Field	
Sensitivity	0.1310	'#2	
Sample Limit 1	-		
Sample Limit 2	-		

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.40	0	
1	921	P	Primer	2.36	1	
2	921	D	Drift	2.29	1	
3	921	C	2.22	2.23	1	
4	901	C	1.66	1.68	1	
5	902	C	1.11	1.06	1	
6	903	C	0.55	0.50	1	
7	904	C	0.22	0.32	1	
8	905	C	Blank	-0.04	1	
9	0	B	Baseline	-0.40	1	
10	921	H1	High	2.45	1	
11	904	L1	Low	0.48	1	
12	904	L1	Low	0.48	1	
13	905	S	Blank	0.08	1	
14	906	K	Blank Spike	1.00	1	
15	907	QC1	CRM [1.5 mg/L]	1.48	1	
16	0	B	Baseline	-0.40	1	
17	1	S	608	0.35	1	
18	2	S	673	0.16	1	
19	3	S	361	0.15	1	
20	4	S	660	0.28	1	
21	5	S	452	-0.17	1	
22	6	S	426	-0.17	1	
23	7	S	621	0.01	1	
24	8	S	647	-0.19	1	
25	9	S	348	0.02	1	
26	10	K	348 spiked	1.11	1	
27	907	QC1	CRM [1.5 mg/L]	1.40	1	
28	0	B	Baseline	-0.40	1	
29	11	S	686	-0.23	1	
30	12	S	530	0.06	1	
31	13	S	491	0.12	1	
32	14	S	634	0.10	1	
33	15	S	374	0.00	1	
34	16	K	374 spiked	-0.16	1	
35	17	S	543	-0.09	1	
36	18	K	543 spiked	1.06	1	
37	907	QC1	CRM [1.5 mg/L]	1.36	1	
38	0	B	Baseline	-0.40	1	
39	19	S	36	-0.12	1	
40	20	S	10	-0.05	1	
41	21	S	192	0.08	1	
42	22	S	270	-0.09	1	

43	23	S	49	-0.05	1
44	24	S	218	-0.04	1
45	25	S	114	-0.22	1
46	26	S	166	-0.02	1
47	27	S	205	-0.10	1
48	28	K	205 spiked	1.03	1
49	907	QC1	CRM [1.5 mg/L]	1.45	1
50	0	B	Baseline	-0.40	1
51	29	S	62	-0.17	1
52	30	S	127	-0.06	1
53	31	S	153	0.04	1
54	32	S	23	0.32	1
55	33	S	101	-0.03	1
56	34	S	88	-0.22	1
57	35	S	140	0.02	1
58	36	S	75	-0.08	1
59	37	S	257	-0.32	1
60	38	S	179	0.01	1
61	39	K	179 spiked	1.02	1
62	907	QC1	CRM [1.5 mg/L]	1.39	1
63	0	B	Baseline	-0.40	1
64	905	S	Blank	-0.19	1
65	905	S	Blank	-0.15	1
66	905	S	Blank	-0.13	1
67	905	S	Blank	-0.15	1
68	905	S	Blank	-0.17	1
69	905	S	Blank	-0.19	1
70	905	S	Blank	-0.20	1
71	905	S	Blank	-0.21	1
72	905	S	Blank	-0.20	1
73	905	S	Blank	-0.22	1
74	906	K	Blank spiked	0.83	1
75	907	QC1	CRM [1.5 mg/L]	1.30	1
76	0	B	Baseline	-0.40	1
77	921	D	Drift	2.29	1
78	0	B	Final Base	-0.40	1

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	done
%	3.67

Spike Recovery


Blank Spike Added	1.00
Blank Spike Found	0.92
Recovery %	92.14
348 spiked Added	1.00
348 spiked Found	1.08
Recovery %	108.38
374 spiked Added	1.00
374 spiked Found	-0.16
Recovery %	-15.63
543 spiked Added	1.00
543 spiked Found	1.15
Recovery %	114.70
205 spiked Added	1.00
205 spiked Found	1.14
Recovery %	113.44
179 spiked Added	1.00
179 spiked Found	1.01
Recovery %	100.88
Blank spiked Added	1.00
Blank spiked Found	1.05
Recovery %	104.59

Calibration

Channel 1 TKN

Type	Observed	Calculated	Target	% Error
1C	57439	2.23	2.22	0.57
2C	46203	1.68	1.66	0.97
3C	33849	1.06	1.11	-4.14
4C	22533	0.50	0.55	-8.47
5C	18820	0.32	0.22	45.21
6C	11653	-0.04	0.00	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9979
 Carryover (%): 3.6697
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 1.237165E+004
 b = 2.018526E+004

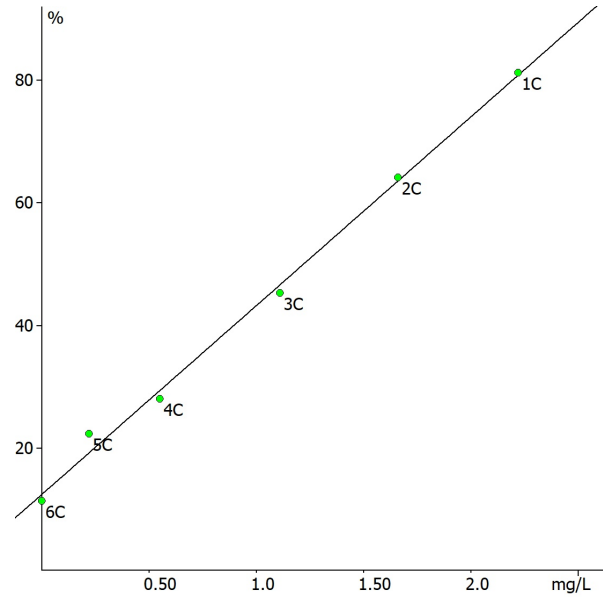

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero