

General

Name of Run	250306_Q8_NP2.run	Name of Analysis	NOx_P04.ANL
Date of Report	3/6/2025	System	AA500
Run Start	3/6/2025 11:08:20 AM	Run Stop	12:26:14 PM
Operator Run/Rec.	CLMR&CVMD	Software Version	AACE 8.02 SP2
Comment			

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	1.0000	0.9985	
AUFS/Gain	0.40/25	0.10/100	Data
Offset	235	5013	Field
Sensitivity	0.3349	0.0830	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.004	0.002	0
1	921	P	Primer	0.203	0.202	1
2	921	D	Drift	0.199	0.204	1
3	921	C	0.200	0.200	0.204	1
4	901	C	0.100	0.101	0.094	1
5	902	C	0.050	0.051	0.045	1
6	903	C	0.025	0.025	0.025	1
7	904	C	0.010	0.010	0.011	1
8	905	C	0.005	0.005	0.006	1
9	906	C	Blank	-0.001	0.004	1
10	0	B	Baseline	-0.004	0.002	1
11	921	H1	High	0.201	0.203	1
12	905	L1	Low	0.006	0.007	1
13	905	L1	Low	0.006	0.007	1
14	906	S	Blank	-0.001	0.004	1
15	907	K	Blank spiked	0.151	0.154	1
16	908	QC1	CRM 2/3 [0.15/0.03]	0.151	0.033	1
17	0	B	Baseline	-0.004	0.002	1
18	1	S	152	0.008	0.012	1
19	2	S	48	0.007	0.016	1
20	3	S	61	0.004	0.014	1
21	4	S	204	0.009	0.003	1
22	5	S	113	0.016	0.003	1
23	6	S	87	0.003	0.011	1
24	7	K	87 spiked	0.161	0.172	1
25	908	QC1	CRM 2/3 [0.15/0.03]	0.151	0.032	1
26	0	B	Baseline	-0.004	0.002	1
27	8	S	22	0.003	0.013	1
28	9	S	9	0.008	0.011	1
29	10	S	74	0.003	0.020	1
30	11	S	100	0.004	0.012	1
31	12	S	139	0.004	0.011	1
32	13	S	35	0.004	0.011	1
33	14	K	35 spiked	0.163	0.152	1
34	908	QC1	CRM 2/3 [0.15/0.03]	0.150	0.030	1
35	0	B	Baseline	-0.004	0.002	1
36	906	S	Blank	-0.001	0.003	1
37	906	S	Blank	-0.001	0.003	1
38	906	S	Blank	-0.001	0.003	1
39	906	S	Blank	-0.001	0.003	1
40	906	S	Blank	-0.001	0.003	1
41	906	S	Blank	-0.001	0.003	1
42	906	S	Blank	-0.001	0.003	1

43	906	S	Blank	-0.001	0.003	1
44	906	S	Blank	-0.001	0.003	1
45	906	S	Blank	-0.001	0.003	1
46	907	K	Blank spiked	0.148	0.152	1
47	908	QC1	CRM 2/3 [0.15/0.03]	0.150	0.033	1
48	0	B	Baseline	-0.004	0.002	1
49	921	D	Drift	0.199	0.204	1
50	0	B	Final Base	-0.004	0.002	1

Corrections

Channel	1	3
Baseline	done	done
Drift	done	done
Carryover	done	done
%	0.87	0.45

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.152	0.151
Recovery %	100.997	100.529
87 spiked Added	0.150	0.150
87 spiked Found	0.157	0.161
Recovery %	104.831	107.695
35 spiked Added	0.150	0.150
35 spiked Found	0.159	0.141
Recovery %	106.125	94.493
Blank spiked Added	0.150	0.150
Blank spiked Found	0.149	0.149
Recovery %	99.301	99.521

Calibration

Channel 1	NOx
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	57190	0.200	0.200	-0.22
2C	30959	0.101	0.100	0.56
3C	17744	0.051	0.050	1.36
4C	10950	0.025	0.025	0.15
5C	6852	0.010	0.010	-4.28
6C	5738	0.005	0.005	7.34
7C	4111	-0.001	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 1.0000
 Carryover (%): 0.8670
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 4.315951E+003
 b = 2.649547E+005

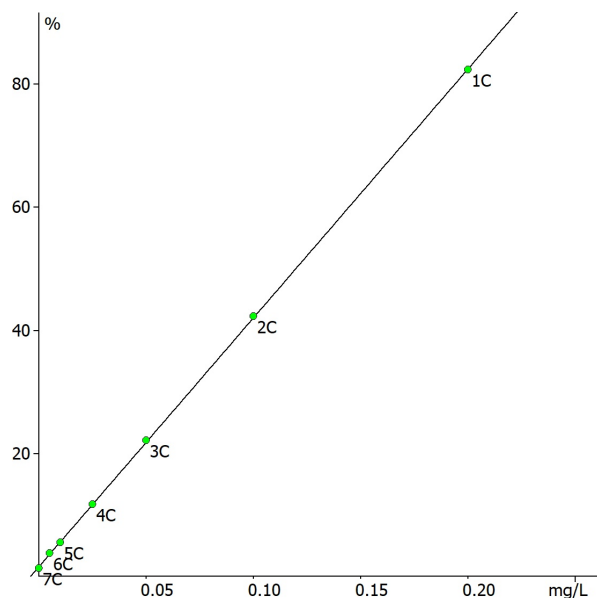

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58033	0.204	0.200	2.04
2C	28141	0.094	0.100	-6.01
3C	14852	0.045	0.050	-9.91
4C	9516	0.025	0.025	1.57
5C	5704	0.011	0.010	13.53
6C	4374	0.006	0.005	29.10
7C	3618	0.004	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9985
 Carryover (%): 0.4466
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 2.621485E+003
 b = 2.715067E+005

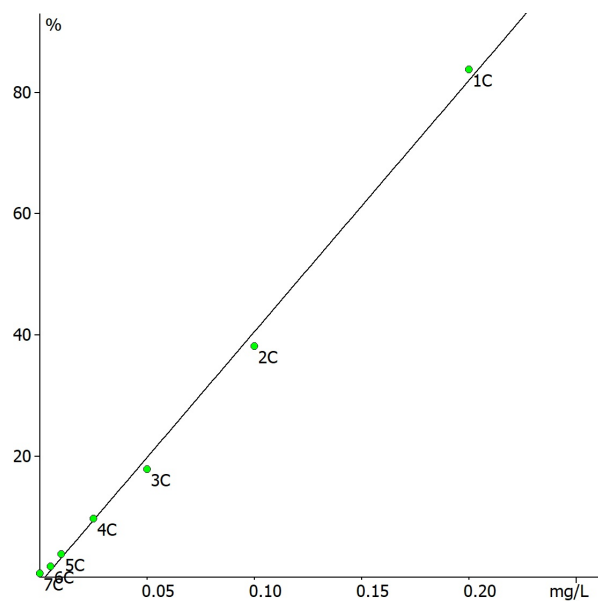

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero