

General

Name of Run	250308_Q8_NP4.run	Name of Analysis	NOx_P04.ANL
Date of Report	4/24/2025	System	AA500
Run Start	3/8/2025 6:23:20 PM	Run Stop	7:39:58 PM
Operator Run/Rec.	CVMD&CLMR	Software Version	AACE 8.02 SP2
Comment			

250308_Q8_NP3.run

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	1.0000	0.9963	
AUFS/Gain	0.42/24	0.10/101	Data
Offset	230	4984	Field
Sensitivity	0.3459	0.0824	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.004	0.003	0
1	921	P	Primer	0.197	0.201	1
2	921	D	Drift	0.198	0.202	1
3	906	C	Blank	0.002	0.006	1
4	0	B	Baseline	-0.004	0.003	1
5	921	H1	High	0.198	0.201	1
6	905	L1	Low	0.004	0.009	1
7	905	L1	Low	0.004	0.009	1
8	906	S	Blank	0.002	0.005	1
9	907	K	Blank spiked	0.149	0.150	1
10	908	QC1	CRM 2/3 [0.15/0.03]	0.149	0.034	1
11	0	B	Baseline	-0.004	0.003	1
12	1	S	126	0.002	0.008	1
13	2	S	516	0.002	0.006	1
14	3	S	334	0.001	0.009	1
15	4	S	178	-0.002	0.007	1
16	5	S	399	0.005	0.006	1
17	6	S	555	0.002	0.005	1
18	7	S	256	0.000	0.006	1
19	8	S	464	-0.001	0.005	1
20	9	S	568	0.002	0.005	1
21	10	S	594	0.001	0.015	1
22	11	K	594 spiked	0.158	0.163	1
23	908	QC1	CRM 2/3 [0.15/0.03]	0.149	0.033	1
24	0	B	Baseline	-0.004	0.003	1
25	12	S	438	0.001	0.005	1
26	13	S	503	-0.001	0.005	1
27	14	S	412	0.003	0.010	1
28	15	S	165	0.003	0.004	1
29	16	S	581	0.003	0.005	1
30	17	S	191	0.003	0.006	1
31	18	S	321	0.000	0.005	1
32	19	K	321 spiked	0.149	0.156	1
33	908	QC1	CRM 2/3 [0.15/0.03]	0.148	0.033	1
34	0	B	Baseline	-0.004	0.003	1
35	906	S	Blank	0.001	0.004	1
36	906	S	Blank	0.001	0.004	1
37	906	S	Blank	0.001	0.004	1
38	906	S	Blank	0.001	0.004	1
39	906	S	Blank	0.001	0.004	1
40	906	S	Blank	0.000	0.005	1
41	906	S	Blank	0.000	0.005	1

42	906	S	Blank	0.000	0.004	1
43	906	S	Blank	0.000	0.005	1
44	906	S	Blank	0.000	0.004	1
45	907	K	Blank spiked	0.148	0.147	1
46	908	QC1	CRM 2/3 [0.15/0.03]	0.149	0.034	1
47	0	B	Baseline	-0.004	0.003	1
48	921	D	Drift	0.198	0.202	1
49	0	B	Final Base	-0.004	0.003	1

Corrections

Channel	1	3
Baseline	done	done
Drift	done	done
Carryover	done	done
%	0.26	0.40

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.146	0.145
Recovery %	97.579	96.515
594 spiked Added	0.150	0.150
594 spiked Found	0.157	0.149
Recovery %	104.708	99.327
321 spiked Added	0.150	0.150
321 spiked Found	0.149	0.150
Recovery %	99.490	100.434
Blank spiked Added	0.150	0.150
Blank spiked Found	0.148	0.143
Recovery %	98.842	95.537

Calibration

Channel 1	NOx
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	58230	0.199	0.200	-0.29
2C	31593	0.101	0.100	0.75
3C	18059	0.051	0.050	1.24
4C	11269	0.025	0.025	1.88
5C	7038	0.010	0.010	-2.04
6C	5714	0.005	0.005	-2.16
7C	4136	-0.001	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 1.0000
 Carryover (%): 0.2642
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 4.393365E+003
 b = 2.699627E+005

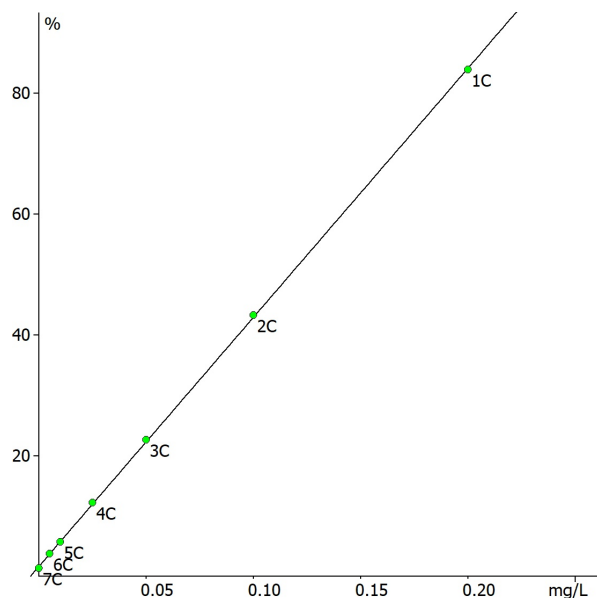

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58616	0.204	0.200	2.12
2C	29181	0.097	0.100	-2.73
3C	12611	0.037	0.050	-25.88
4C	9709	0.027	0.025	6.05
5C	5555	0.011	0.010	14.18
6C	4859	0.009	0.005	77.78
7C	3680	0.005	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9963
 Carryover (%): 0.3974
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 2.412925E+003
 b = 2.751808E+005

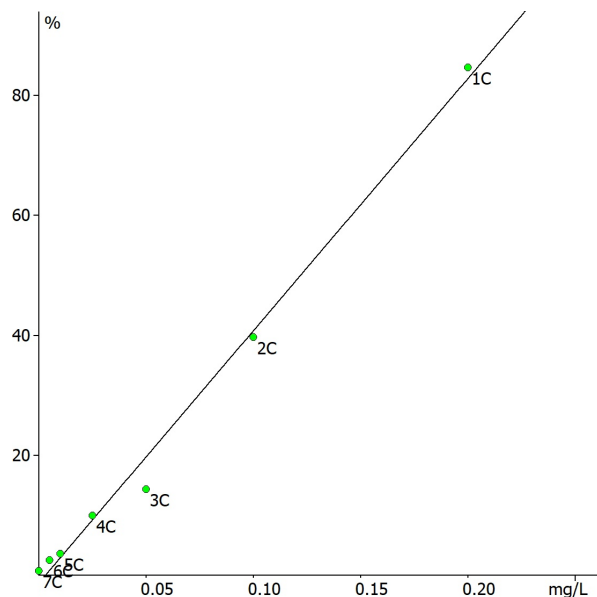

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero