

## General

|                   |                      |                  |               |
|-------------------|----------------------|------------------|---------------|
| Name of Run       | 250312_Q8_TKN1.run   | Name of Analysis | TKN.ANL       |
| Date of Report    | 4/24/2025            | System           | AA500         |
| Run Start         | 3/12/2025 1:58:29 PM | Run Stop         | 3:56:02 PM    |
| Operator Run/Rec. | CLMR&CVMD            | Software Version | AACE 8.02 SP2 |
| Comment           | 2nd Batch            |                  |               |

241206\_Q7\_TKN1.run

## Parameters

|                |         |          |  |
|----------------|---------|----------|--|
| Channel        | 1       | Dilution |  |
| Method         | TKN     |          |  |
| Unit           | mg/L    |          |  |
| Calib. Fit     | Linear  |          |  |
| Corr.Coeff.(r) | 0.9974  |          |  |
| AUFS/Gain      | 0.15/66 | Data     |  |
| Offset         | -2330   | Field    |  |
| Sensitivity    | 0.1417  | '#2      |  |
| Sample Limit 1 | -       |          |  |
| Sample Limit 2 | -       |          |  |

## Results

| Pk | Cup | Type | Sample ID      | Value | Data | Value |
|----|-----|------|----------------|-------|------|-------|
| 0  | 0   | B    | Initial Base   | -0.43 | 0    |       |
| 1  | 921 | P    | Primer         | 2.28  | 1    |       |
| 2  | 921 | D    | Drift          | 2.36  | 1    |       |
| 3  | 921 | C    | 2.22           | 2.27  | 1    |       |
| 4  | 901 | C    | 1.66           | 1.66  | 1    |       |
| 5  | 902 | C    | 1.11           | 1.06  | 1    |       |
| 6  | 903 | C    | 0.55           | 0.45  | 1    |       |
| 7  | 904 | C    | 0.22           | 0.25  | 1    |       |
| 8  | 905 | C    | Blank          | 0.07  | 1    |       |
| 9  | 0   | B    | Baseline       | -0.43 | 1    |       |
| 10 | 921 | H1   | High           | 2.11  | 1    |       |
| 11 | 904 | L1   | Low            | 0.37  | 1    |       |
| 12 | 904 | L1   | Low            | 0.40  | 1    |       |
| 13 | 905 | S    | Blank          | 0.16  | 1    |       |
| 14 | 906 | K    | Blank Spike    | 1.09  | 1    |       |
| 15 | 907 | QC1  | CRM [1.5 mg/L] | 1.63  | 1    |       |
| 16 | 0   | B    | Baseline       | -0.43 | 1    |       |
| 17 | 1   | S    | 335            | 0.11  | 1    |       |
| 18 | 2   | S    | 400            | 0.19  | 1    |       |
| 19 | 3   | S    | 101            | 0.21  | 1    |       |
| 20 | 4   | S    | 62             | 0.15  | 1    |       |
| 21 | 5   | S    | 127            | 0.21  | 1    |       |
| 22 | 6   | S    | 517            | 0.22  | 1    |       |
| 23 | 7   | S    | 413            | 0.25  | 1    |       |
| 24 | 8   | S    | 595            | 0.08  | 1    |       |
| 25 | 9   | S    | 140            | 0.07  | 1    |       |
| 26 | 10  | S    | 10             | 0.02  | 1    |       |
| 27 | 11  | K    | 10 spiked      | 1.13  | 1    |       |
| 28 | 907 | QC1  | CRM [1.5 mg/L] | 1.61  | 1    |       |
| 29 | 0   | B    | Baseline       | -0.43 | 1    |       |
| 30 | 12  | S    | 166            | -0.01 | 1    |       |
| 31 | 13  | S    | 257            | 0.05  | 1    |       |
| 32 | 14  | S    | 23             | 0.24  | 1    |       |
| 33 | 15  | S    | 179            | 0.10  | 1    |       |
| 34 | 16  | S    | 36             | 0.12  | 1    |       |
| 35 | 17  | S    | 205            | 0.13  | 1    |       |
| 36 | 18  | S    | 465            | -0.01 | 1    |       |
| 37 | 19  | S    | 566            | -0.01 | 1    |       |
| 38 | 20  | S    | 114            | 0.01  | 1    |       |
| 39 | 21  | S    | 153            | 0.03  | 1    |       |
| 40 | 22  | K    | 153 spiked     | 1.07  | 1    |       |
| 41 | 907 | QC1  | CRM [1.5 mg/L] | 1.60  | 1    |       |
| 42 | 0   | B    | Baseline       | -0.43 | 1    |       |

|    |     |     |                |       |   |
|----|-----|-----|----------------|-------|---|
| 43 | 23  | S   | 439            | 0.01  | 1 |
| 44 | 24  | S   | 192            | -0.05 | 1 |
| 45 | 25  | S   | 582            | 0.06  | 1 |
| 46 | 26  | S   | 49             | -0.17 | 1 |
| 47 | 27  | S   | 88             | 0.06  | 1 |
| 48 | 28  | S   | 569            | 0.07  | 1 |
| 49 | 29  | S   | 504            | 0.04  | 1 |
| 50 | 30  | S   | 322            | 0.17  | 1 |
| 51 | 31  | S   | 75             | 0.02  | 1 |
| 52 | 32  | K   | 75 spiked      | 1.11  | 1 |
| 53 | 907 | QC1 | CRM [1.5 mg/L] | 1.62  | 1 |
| 54 | 0   | B   | Baseline       | -0.43 | 1 |
| 55 | 905 | S   | Blank          | -0.06 | 1 |
| 56 | 905 | S   | Blank          | -0.08 | 1 |
| 57 | 905 | S   | Blank          | -0.05 | 1 |
| 58 | 905 | S   | Blank          | -0.03 | 1 |
| 59 | 905 | S   | Blank          | -0.03 | 1 |
| 60 | 905 | S   | Blank          | -0.05 | 1 |
| 61 | 905 | S   | Blank          | -0.06 | 1 |
| 62 | 905 | S   | Blank          | -0.07 | 1 |
| 63 | 905 | S   | Blank          | -0.09 | 1 |
| 64 | 905 | S   | Blank          | -0.12 | 1 |
| 65 | 906 | K   | Blank spiked   | 0.96  | 1 |
| 66 | 907 | QC1 | CRM [1.5 mg/L] | 1.62  | 1 |
| 67 | 0   | B   | Baseline       | -0.43 | 1 |
| 68 | 921 | D   | Drift          | 2.36  | 1 |
| 69 | 0   | B   | Final Base     | -0.43 | 1 |

## Corrections

|           |       |
|-----------|-------|
| Channel   | 1     |
| Baseline  | done  |
| Drift     | done  |
| Carryover | none  |
| %         | -1.34 |

## Spike Recovery


|                    |        |
|--------------------|--------|
| Blank Spike Added  | 1.00   |
| Blank Spike Found  | 0.93   |
| Recovery %         | 93.42  |
| 10 spiked Added    | 1.00   |
| 10 spiked Found    | 1.11   |
| Recovery %         | 111.06 |
| 153 spiked Added   | 1.00   |
| 153 spiked Found   | 1.04   |
| Recovery %         | 104.28 |
| 75 spiked Added    | 1.00   |
| 75 spiked Found    | 1.09   |
| Recovery %         | 108.99 |
| Blank spiked Added | 1.00   |
| Blank spiked Found | 1.09   |
| Recovery %         | 108.57 |

## Calibration

Channel 1 TKN

| Type | Observed | Calculated | Target | % Error |
|------|----------|------------|--------|---------|
| 1C   | 63326    | 2.27       | 2.22   | 2.03    |
| 2C   | 49577    | 1.66       | 1.66   | -0.07   |
| 3C   | 36078    | 1.06       | 1.11   | -4.18   |
| 4C   | 22191    | 0.45       | 0.55   | -17.96  |
| 5C   | 17715    | 0.25       | 0.22   | 15.38   |
| 6C   | 13485    | 0.07       | 0.00   | -       |

Calibration Fit: Linear  
 Correlation Coefficient: 0.9974  
 Carryover (%): no  
 Calibration Equation:  $y = bx + a$   
 $y =$  Peak Height Digit. Units  
 $x =$  Concentration Unit mg/L  
 $a =$  1.195866E+004  
 $b =$  2.267709E+004

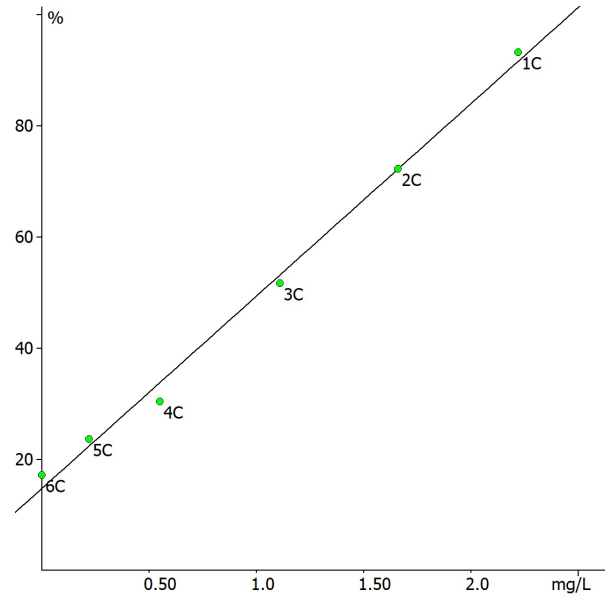

## Legend

\* ... Sample offscale  
 + ... Result higher than sample limit  
 - ... Result lower than sample limit  
 P ... Standard passed  
 F ... Standard failed  
 N ... Value not calculated or not used  
 R ... Resample after offscale  
 M ... Peak marker moved manually  
 d ... Pre-Diluted sample  
 D ... Post-Diluted sample  
 H ... Dual calibration used  
 Z ... Report negative values as zero