

General

Name of Run	250402_Q8_NP5.run	Name of Analysis	NOx_P04.ANL
Date of Report	4/2/2025	System	AA500
Run Start	4/2/2025 2:29:23 PM	Run Stop	3:59:55 PM
Operator Run/Rec. CLMR		Software Version	AACE 8.02 SP2
Comment			

250308_Q8_NP3.run

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9998	1.0000	
AUFS/Gain	0.26/38	0.10/98	Data
Offset	246	5017	Field
Sensitivity	0.2499	0.0851	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.004	-0.001	0
1	921	P	Primer	0.178	0.198	1
2	921	D	Drift	0.197 *	0.199	1
3	921	C	0.200	0.199 *	0.200	1
4	901	C	0.100	0.102	0.100	1
5	902	C	0.050	0.051	0.050	1
6	903	C	0.025	0.025	0.025	1
7	904	C	0.010	0.012	0.010	1
8	905	C	0.005	0.004	0.006	1
9	906	C	Blank	-0.002	0.000	1
10	0	B	Baseline	-0.004	-0.001	1
11	921	H1	High	0.211 *	0.199	1
12	905	L1	Low	0.005	0.006	1
13	905	L1	Low	0.004	0.006	1
14	906	S	Blank	-0.002	0.000	1
15	907	K	Blank spiked	0.166	0.147	1
16	908	QC1	CRM 2/3 [0.15/0.03]	0.162	0.029	1
17	0	B	Baseline	-0.004	-0.001	1
18	1	S	646	0.000	0.002	1
19	2	S	295	-0.001	0.006	1
20	3	S	607	0.000	0.001	1
21	4	S	672	0.000	0.002	1
22	5	S	542	0.000	0.003	1
23	6	S	620	0.002	0.002	1
24	7	S	295	-0.001	0.009	1
25	8	S	308	-0.001	0.003	1
26	9	S	373	-0.002	0.011	1
27	10	S	256	0.001	-0.001	1
28	11	S	269	0.000	0.009	1
29	12	K	269 spiked	0.167	0.157	1
30	908	QC1	CRM 2/3 [0.15/0.03]	0.164	0.027	1
31	0	B	Baseline	-0.004	-0.001	1
32	13	S	464	-0.002	0.000	1
33	14	S	490	0.004	0.005	1
34	15	S	347	0.000	0.001	1
35	16	S	360	0.003	0.003	1
36	17	S	529	0.000	0.002	1
37	18	S	659	0.002	0.002	1
38	19	S	503	0.001	-0.003	1
39	20	S	282	0.006	0.002	1
40	21	S	685	0.000	0.001	1
41	22	S	386	0.001	0.103	1
42	23	S	178	0.000	0.457 *	1

43	24	K	178 spiked	0.157	0.591 *	1
44	908	QC1	CRM 2/3 [0.15/0.03]	0.161	0.027	1
45	0	B	Baseline	-0.004	-0.001	1
46	906	S	Blank	0.000	-0.001	1
47	906	S	Blank	0.000	0.000	1
48	906	S	Blank	0.000	-0.001	1
49	906	S	Blank	0.000	-0.001	1
50	906	S	Blank	0.001	0.000	1
51	906	S	Blank	0.001	0.000	1
52	906	S	Blank	0.000	0.000	1
53	906	S	Blank	0.000	0.000	1
54	906	S	Blank	-0.001	-0.001	1
55	906	S	Blank	-0.001	-0.001	1
56	907	K	Blank spiked	0.168	0.144	1
57	908	QC1	CRM 2/3 [0.15/0.03]	0.166	0.028	1
58	0	B	Baseline	-0.004	-0.001	1
59	921	D	Drift	0.218 *	0.199	1
60	0	B	Final Base	-0.004	-0.001	1

Corrections

Channel	1	3
Baseline	done	done
Drift	none	done
Carryover	none	done
%	0.00	0.49

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.168	0.148
Recovery %	111.744	98.434
269 spiked Added	0.150	0.150
269 spiked Found	0.167	0.147
Recovery %	111.264	98.310
178 spiked Added	0.150	0.136
178 spiked Found	0.158	0.134
Recovery %	105.279	98.547
Blank spiked Added	0.150	0.150
Blank spiked Found	0.169	0.145
Recovery %	112.804	96.891

Calibration

Channel 1	NOx
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	72787	0.199	0.200	-0.60
2C	39433	0.102	0.100	1.64
3C	22160	0.051	0.050	2.63
4C	13100	0.025	0.025	-0.32
5C	8600	0.012	0.010	18.10
6C	5816	0.004	0.005	-26.02
7C	3797	-0.002	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 0.9998
 Carryover (%): 0.0000
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 4.546354E+003
 b = 3.432526E+005

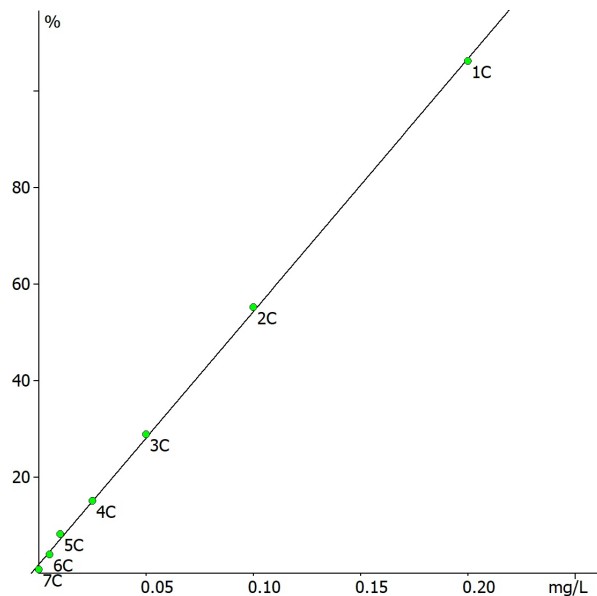

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58432	0.200	0.200	0.15
2C	30844	0.100	0.100	-0.37
3C	17146	0.050	0.050	-0.72
4C	10274	0.025	0.025	-1.74
5C	6244	0.010	0.010	-1.42
6C	5287	0.006	0.005	27.31
7C	3444	0.000	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 1.0000
 Carryover (%): 0.4913
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 3.542665E+003
 b = 2.740317E+005

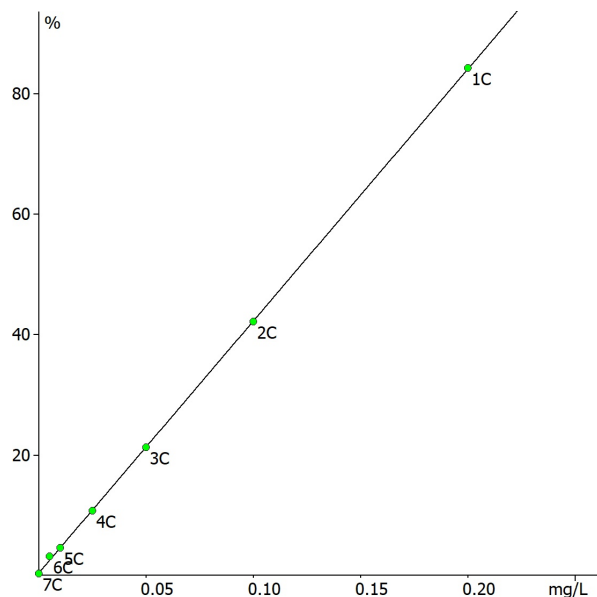

Legend

* ... Sample offscale
 + ... Result higher than sample limit
 - ... Result lower than sample limit
 P ... Standard passed
 F ... Standard failed
 N ... Value not calculated or not used
 R ... Resample after offscale
 M ... Peak marker moved manually
 d ... Pre-Diluted sample
 D ... Post-Diluted sample
 H ... Dual calibration used
 Z ... Report negative values as zero