

General

Name of Run	250407_Q8_TKN3_4CNull.run	Name of Analysis	TKN.ANL
Date of Report	4/24/2025	System	AA500
Run Start	4/7/2025 2:17:08 PM	Run Stop	3:57:27 PM
Operator Run/Rec.	CLMR / CLMR	Software Version	AACE 8.02 SP2
Comment	Recalculate from Run250407_Q8_TKN3.run		

241206_Q7_TKN1.run

Parameters

Channel	1	
Method	TKN	Dilution
Unit	mg/L	
Calib. Fit	Linear	
Corr.Coeff.(r)	0.9982	
AUFS/Gain	0.18/57	Data
Offset	-1436	Field
Sensitivity	0.1495	'#2
Sample Limit 1	-	
Sample Limit 2	-	

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.32	0	
1	921	P	Primer	2.29	1	
2	921	D	Drift	2.30	1	
3	921	C	2.22	2.27	1	
4	901	C	1.66	1.63	1	
5	902	C	1.11	1.05	1	
6	0	X		0.84 N	1	
7	904	C	0.22	0.29	1	
8	905	C	Blank	-0.02	1	
9	0	B	Baseline	-0.32	1	
10	921	H1	High	2.07	1	
11	904	L1	Low	0.27	1	
12	904	L1	Low	0.29	1	
13	905	S	Blank	0.00	1	
14	906	K	Blank Spike	0.91	1	
15	907	QC1	CRM [1.5 mg/L]	1.35	1	
16	0	B	Baseline	-0.32	1	
17	1	S	218	0.38	1	
18	2	S	673	0.19	1	
19	3	S	387	0.19	1	
20	4	S	309	0.07	1	
21	5	S	660	0.04	1	
22	6	S	374	0.08	1	
23	7	S	621	0.02	1	
24	8	S	361	0.01	1	
25	9	S	348	0.04	1	
26	10	K	348 spiked	0.98	1	
27	907	QC1	CRM [1.5 mg/L]	1.31	1	
28	0	B	Baseline	-0.32	1	
29	11	S	686	0.00	1	
30	12	S	283	0.04	1	
31	13	S	543	0.07	1	
32	14	S	608	0.59	1	
33	15	S	491	0.02	1	
34	16	S	270	0.02	1	
35	17	S	530	0.01	1	
36	18	S	647	0.03	1	
37	19	S	296	0.07	1	
38	20	K	296 spiked	1.07	1	
39	907	QC1	CRM [1.5 mg/L]	1.30	1	
40	0	B	Baseline	-0.32	1	
41	905	S	Blank	-0.09	1	
42	905	S	Blank	-0.11	1	

43	905	S	Blank	-0.12	1
44	905	S	Blank	-0.11	1
45	905	S	Blank	-0.09	1
46	905	S	Blank	-0.09	1
47	905	S	Blank	-0.12	1
48	905	S	Blank	-0.07	1
49	905	S	Blank	-0.07	1
50	905	S	Blank	-0.07	1
51	906	K	Blank spiked	0.92	1
52	907	QC1	CRM [1.5 mg/L]	1.38	1
53	0	B	Baseline	-0.32	1
54	921	D	Drift	2.30	1
55	0	B	Final Base	-0.32	1

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	none
%	-1.34

Spike Recovery


Blank Spike Added	1.00
Blank Spike Found	0.92
Recovery %	91.49
348 spiked Added	1.00
348 spiked Found	0.94
Recovery %	94.36
296 spiked Added	1.00
296 spiked Found	1.00
Recovery %	99.69
Blank spiked Added	1.00
Blank spiked Found	0.99
Recovery %	98.81

Calibration

Channel 1	TKN
-----------	-----

Type	Observed	Calculated	Target	% Error
1C	60266	2.27	2.22	2.20
2C	46519	1.63	1.66	-1.97
3C	34085	1.05	1.11	-5.67
5C	17886	0.29	0.22	32.34
6C	11124	-0.02	0.00	-

Calibration Fit: Linear
 Correlation Coefficient: 0.9982
 Carryover (%): no
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 1.164690E+004
 b = 2.142908E+004

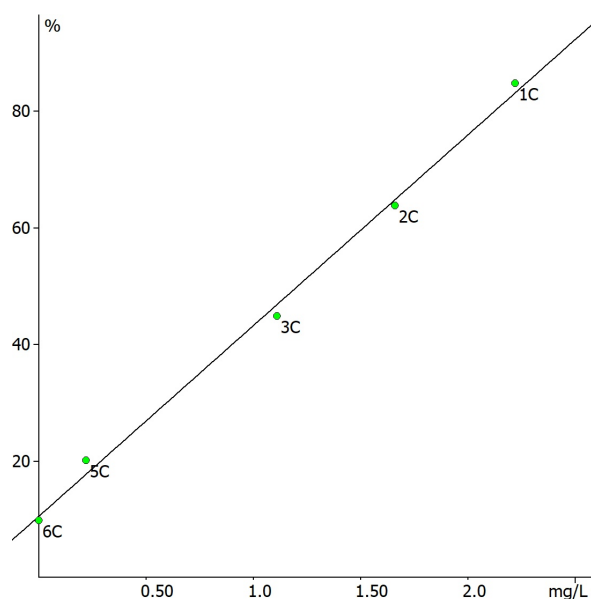

Legend

* ... Sample offscale
+ ... Result higher than sample limit
- ... Result lower than sample limit
P ... Standard passed
F ... Standard failed
N ... Value not calculated or not used
R ... Resample after offscale
M ... Peak marker moved manually
d ... Pre-Diluted sample
D ... Post-Diluted sample
H ... Dual calibration used
Z ... Report negative values as zero