

General

Name of Run	250409_Q8_NP6.run	Name of Analysis	NOx_P04.ANL
Date of Report	4/24/2025	System	AA500
Run Start	4/9/2025 2:01:26 PM	Run Stop	3:04:30 PM
Operator Run/Rec.	CLMR&CVMD	Software Version	AACE 8.02 SP2
Comment			

250308_Q8_NP3.run

Parameters

Channel	1	3	
Method	NOx	P04	Dilution
Unit	mg/L	mg/L	
Calib. Fit	Linear	Linear	
Corr.Coeff.(r)	0.9994	1.0000	
AUFS/Gain	0.23/43	0.10/99	Data
Offset	182	5029	Field
Sensitivity	0.2205	0.0852	'#2
Sample Limit 1	-	-	
Sample Limit 2	-	-	

Results

Pk	Cup	Type	Sample ID	Value	Value	Data
0	0	B	Initial Base	-0.005	-0.002	0
1	921	P	Primer	0.161	0.199	1
2	921	D	Drift	0.193 *	0.199	1
3	921	C	0.200	0.198 *	0.200	1
4	901	C	0.100	0.105	0.100	1
5	902	C	0.050	0.051	0.050	1
6	903	C	0.025	0.024	0.025	1
7	904	C	0.010	0.010	0.010	1
8	905	C	0.005	0.006	0.005	1
9	906	C	Blank	-0.003	0.000	1
10	0	B	Baseline	-0.005	-0.002	1
11	921	H1	High	0.217 *	0.199	1
12	905	L1	Low	0.007	0.005	1
13	905	L1	Low	0.006	0.005	1
14	906	S	Blank	-0.003	-0.001	1
15	907	K	Blank spiked	0.161 *	0.148	1
16	908	QC1	CRM 2/3 [0.15/0.03]	0.167 *	0.029	1
17	0	B	Baseline	-0.005	-0.002	1
18	1	S	711	-0.003	-0.001	1
19	2	S	724	-0.001	-0.001	1
20	3	S	698	0.002	0.000	1
21	4	S	451	-0.002	-0.001	1
22	5	K	451 spiked	0.161 *	0.149	1
23	908	QC1	CRM 2/3 [0.15/0.03]	0.167 *	0.029	1
24	0	B	Baseline	-0.005	-0.002	1
25	906	S	Blank	-0.003	-0.001	1
26	906	S	Blank	-0.003	-0.001	1
27	906	S	Blank	-0.003	-0.001	1
28	906	S	Blank	-0.003	-0.001	1
29	906	S	Blank	-0.003	-0.001	1
30	906	S	Blank	-0.003	-0.001	1
31	906	S	Blank	-0.003	0.000	1
32	906	S	Blank	-0.003	-0.001	1
33	906	S	Blank	-0.003	-0.001	1
34	906	S	Blank	-0.003	-0.001	1
35	907	K	Blank spiked	0.161 *	0.147	1
36	908	QC1	CRM 2/3 [0.15/0.03]	0.167 *	0.029	1
37	0	B	Baseline	-0.005	-0.002	1
38	921	D	Drift	0.221 *	0.199	1
39	0	B	Final Base	-0.005	-0.002	1

Corrections

Channel	1	3
Baseline	done	done
Drift	none	done
Carryover	none	done
%	0.00	0.57

Spike Recovery


Blank spiked Added	0.150	0.150
Blank spiked Found	0.164	0.149
Recovery %	109.292	99.259
451 spiked Added	0.150	0.150
451 spiked Found	0.163	0.149
Recovery %	108.808	99.592
Blank spiked Added	0.150	0.150
Blank spiked Found	0.164	0.148
Recovery %	109.543	98.463

Calibration

Channel 1 NOx

Type	Observed	Calculated	Target	% Error
1C	79247	0.198	0.200	-1.24
2C	44441	0.105	0.100	4.72
3C	24156	0.051	0.050	1.27
4C	14334	0.024	0.025	-2.22
5C	8977	0.010	0.010	1.62
6C	7266	0.006	0.005	11.99
7C	4007	-0.003	0.000	-


Calibration Fit: Linear
 Correlation Coefficient: 0.9994
 Carryover (%): 0.0000
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 5.165867E+003
 b = 3.750434E+005

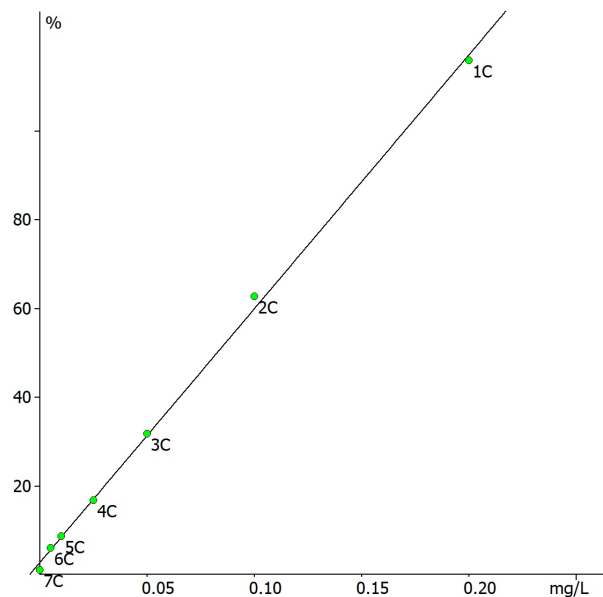

Channel 3 P04

Type	Observed	Calculated	Target	% Error
1C	58750	0.200	0.200	-0.10
2C	31348	0.100	0.100	0.42
3C	17471	0.050	0.050	0.20
4C	10465	0.025	0.025	-1.23
5C	6435	0.010	0.010	0.78
6C	5023	0.005	0.005	-0.85
7C	3642	0.000	0.000	-

Calibration Fit: Linear
 Correlation Coefficient: 1.0000
 Carryover (%): 0.5672
 Calibration Equation: $y = bx + a$
 y = Peak Height Digit. Units
 x = Concentration Unit mg/L
 a = 3.656034E+003
 b = 2.757481E+005

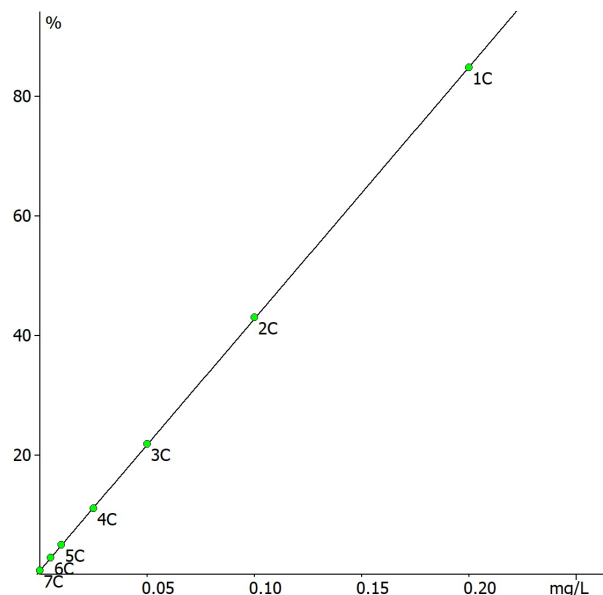

Legend

* ... Sample offscale
+ ... Result higher than sample limit
- ... Result lower than sample limit
P ... Standard passed
F ... Standard failed
N ... Value not calculated or not used
R ... Resample after offscale
M ... Peak marker moved manually
d ... Pre-Diluted sample
D ... Post-Diluted sample
H ... Dual calibration used
Z ... Report negative values as zero