

General

Name of Run	250410_Q8_TKN4_4CNull.run	Name of Analysis	TKN.ANL
Date of Report	4/24/2025	System	AA500
Run Start	4/10/2025 2:44:07 PM	Run Stop	3:54:45 PM
Operator Run/Rec.	CLMR+CVMD / CLMR	Software Version	AACE 8.02 SP2
Comment	Recalculate from Run250410_Q8_TKN4.run		

241206_Q7_TKN1.run

Parameters

Channel	1	
Method	TKN	Dilution
Unit	mg/L	
Calib. Fit	Linear	
Corr.Coeff.(r)	0.9965	
AUFS/Gain	0.21/47	Data
Offset	-2729	Field
Sensitivity	0.1811	'#2
Sample Limit 1	-	
Sample Limit 2	-	

Results

Pk	Cup	Type	Sample ID	Value	Data	Value
0	0	B	Initial Base	-0.55	0	
1	921	P	Primer	2.71	1	
2	921	D	Drift	2.50	1	
3	921	C	2.22	2.32	1	
4	901	C	1.66	1.58	1	
5	902	C	1.11	1.03	1	
6	0	X		0.85 N	1	
7	904	C	0.22	0.26	1	
8	905	C	Blank	0.02	1	
9	0	B	Baseline	-0.55	1	
10	921	H1	High	2.34	1	
11	904	L1	Low	0.51	1	
12	904	L1	Low	0.53	1	
13	905	S	Blank	0.13	1	
14	906	K	Blank Spike	1.03	1	
15	907	QC1	CRM [1.5 mg/L]	1.42	1	
16	0	B	Baseline	-0.55	1	
17	1	S	452	0.13	1	
18	2	S	699	0.41	1	
19	3	S	712	0.45	1	
20	4	S	725	0.35	1	
21	5	K	725 spiked	1.25	1	
22	907	QC1	CRM [1.5 mg/L]	1.33	1	
23	0	B	Baseline	-0.55	1	
24	905	S	Blank	-0.06	1	
25	905	S	Blank	0.13	1	
26	905	S	Blank	0.15	1	
27	905	S	Blank	0.07	1	
28	905	S	Blank	0.05	1	
29	905	S	Blank	0.01	1	
30	905	S	Blank	0.08	1	
31	905	S	Blank	0.05	1	
32	905	S	Blank	0.02	1	
33	905	S	Blank	-0.03	1	
34	906	K	Blank spiked	0.93	1	
35	907	QC1	CRM [1.5 mg/L]	1.39	1	
36	0	B	Baseline	-0.55	1	
37	921	D	Drift	2.50	1	
38	0	B	Final Base	-0.55	1	

Corrections

Channel	1
Baseline	done
Drift	done
Carryover	done
%	5.28

Spike Recovery


Blank Spike Added	1.00
Blank Spike Found	0.90
Recovery %	89.85
725 spiked Added	1.00
725 spiked Found	0.91
Recovery %	90.77
Blank spiked Added	1.00
Blank spiked Found	0.95
Recovery %	95.33

Calibration

Channel 1 TKN

Type	Observed	Calculated	Target	% Error
1C	51360	2.32	2.22	4.36
2C	38766	1.58	1.66	-4.78
3C	29282	1.03	1.11	-7.54
5C	16231	0.26	0.22	19.73
6C	12114	0.02	0.00	-

Calibration Fit: Linear
Correlation Coefficient: 0.9965
Carryover (%): 5.2810
Calibration Equation: $y = bx + a$
y = Peak Height Digit. Units
x = Concentration Unit mg/L
a = 1.172478E+004
b = 1.710731E+004

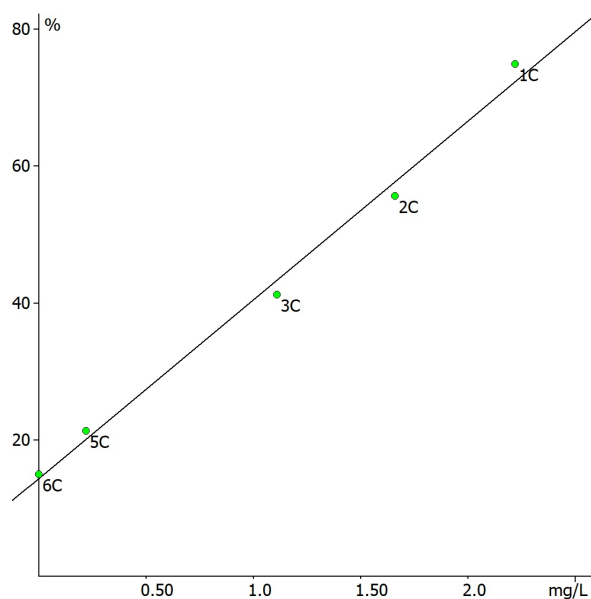

Legend

* ... Sample offscale
+ ... Result higher than sample limit
- ... Result lower than sample limit
P ... Standard passed
F ... Standard failed
N ... Value not calculated or not used
R ... Resample after offscale
M ... Peak marker moved manually
d ... Pre-Diluted sample
D ... Post-Diluted sample
H ... Dual calibration used
Z ... Report negative values as zero